



BMJ Open Performance and safety of PowerPICC catheters and accessories: a prospective observational study

Martin Troubil,¹ Giuseppe Capozzoli,² Baudolino Mussa,³ Melinda Hodne,⁴ Klaus Hoerauf ,^{4,5} Kimberly Alsbrooks ⁴

To cite: Troubil M, Capozzoli G, Mussa B, *et al*. Performance and safety of PowerPICC catheters and accessories: a prospective observational study. *BMJ Open* 2024;**14**:e081288. doi:10.1136/bmjopen-2023-081288

► Prepublication history and additional supplemental material for this paper are available online. To view these files, please visit the journal online (<https://doi.org/10.1136/bmjopen-2023-081288>).

Received 23 October 2023
Accepted 16 August 2024



© Author(s) (or their employer(s)) 2024. Re-use permitted under CC BY-NC. No commercial re-use. See rights and permissions. Published by BMJ.

¹University Hospital Olomouc, Olomouc, Czech Republic

²Bolzano Central Hospital, Bolzano, Italy

³Molinette Hospital, University of Turin, Turin, Italy

⁴Becton, Dickinson and Company, Franklin Lakes, New Jersey, USA

⁵Medical University of Vienna, Vienna, Austria

Correspondence to
Dr Kimberly Alsbrooks;
kim.alsbrooks@bd.com

ABSTRACT

Objective This study aimed to evaluate the safety and performance of PowerPICC catheters in a real-world setting.

Design Prospective, observational, multicentre study.

Setting Nine European countries, involving 14 centres.

Participants General patient population.

Intervention PowerPICC catheter inserted by the clinician as standard of care with routinely collected outcomes followed through device removal or 180 days postinsertion.

Primary and secondary outcomes measures Safety and performance outcomes were assessed for PowerPICC, PowerPICC SOLO 2 and PowerGroshong PICC. The primary safety endpoint was the incidence of symptomatic venous thrombosis (VT), and secondary safety endpoints included phlebitis, extravasation, vessel laceration, vessel perforation local infection, accidental dislodgment and catheter-related bloodstream infection (CRBSI). The primary performance endpoint was the percentage of patients whose PowerPICC device remained in place through the completion of therapy. The secondary performance endpoints included catheter patency, placement success in a single attempt and usability.

Results The enrolled patients (N=451) received either PowerPICC, PowerPICC SOLO 2 or PowerGroshong PICC catheters. Across all devices, 1.6% of patients developed symptomatic VT, and CRBSI occurred in 1.6% of patients. There were no cases of phlebitis or extravasation and only three cases of vein laceration or vein perforation. The catheters showed high success rates in completing therapy (81.8%), maintaining patency (93.9%) and achieving successful placement in a single attempt (90.4%). Clinicians overwhelmingly agreed that both the guidewire and stylet (93.3% and 94.4%, respectively) were easy or very easy to use.

Conclusions This study demonstrates the safety and performance of PowerPICC catheters across diverse settings and patient cohorts in real-world hospital settings across Europe. The findings indicate that these catheters are safe and can be effectively used in the general patient setting and when inserted by a variety of clinicians. The low incidence of complications and high success rates further support the clinical utility of these catheters.

Trial registration number NCT04263649.

INTRODUCTION

Vascular access is a critical component for developing and delivering healthcare plans for hospitalised patients. Approximately 90%

STRENGTHS AND LIMITATIONS OF THIS STUDY

- ⇒ This study provides real-world evidence assessment on the safety and efficacy of PowerPICC catheters.
- ⇒ Observations across 9 European countries and 14 centres revealed high safety and performance, despite potential variations in clinical practice across hospitals.
- ⇒ The study population was diverse, including patients admitted for various health diagnoses such as oncology, infection, surgical, cardiac, gastrointestinal and neurological conditions.
- ⇒ Patients were not randomised and underlying conditions, as well as clinical practices at the specific site, may have impacted the findings.

of patients require vascular access during hospitalisation¹ and millions of vascular access devices (VADs) are placed across Europe every year for the administration of medications, fluids, and blood products as well as for blood sampling, haemodynamic monitoring, and power injection of contrast media.²⁻⁵ These VADs include peripheral intravenous catheters (PIVCs), peripherally inserted central catheters (PICCs), femorally inserted central catheters (FICCs), and centrally inserted central catheters (CICCs).

The selection of a particular VAD should be tailored to the needs of the patient by considering a combination of clinical factors, including the duration of therapy, infusion requirements (eg, need for multiple lumens), and risk of complications.⁶ For patients requiring central access, PICCs are often preferred over CICCs and FICCs due to a reduced risk of complications associated with placement via subclavian, jugular, and femoral veins.⁷ Notably, hospital systems with dedicated PICC placement teams have been shown to confer additional advantages with PICC utilisation by optimising PICC placement; these benefits include decreased patient wait times, improved patient care, and decreased costs associated with PICC

placement.⁸⁹ PICCs may also be employed to decrease the frequency of venipunctures required for the replacement of PIVCs in patients who require therapy for prolonged durations and for administration of infusates with certain osmolarity, pH, or cytotoxicity profiles.

Patients who receive PICCs may also require infusion of contrast media for enhanced CT. Contrast media can enhance the opacity of a target tissue or anatomical structure in order to diagnose and monitor certain diseases, including cancer, arterial stenosis, and pulmonary embolism.^{10 11} However, injection of high-pressure contrast media requires delivery of these agents at well-controlled flow rates with reinforced lumens capable of withstanding the pressures associated with timed contrast infusions.¹² The need for high-pressure contrast media injection-compatible VADs has led to the development of VADs with 'power injection' capabilities. A range of power injectable-compatible PICCs are now available with multipurpose functions that provide vascular access for injection of contrast media in addition to standard uses such as intravenous therapy, blood sampling, parenteral nutrition, and venous pressure measurements.

In particular, the PowerPICC line of devices (BD, Franklin Lakes, New Jersey, USA) was introduced to the European market in 2008 and has been successfully used in hospital systems for over a decade. Therefore, it is important to fully understand the safety and performance of these devices. This is the first prospective study to further evaluate the overall safety and performance of PowerPICC catheters in a real-world setting using an observational, multicentre approach.

METHODS

Study overview and oversight

This prospective, observational, multicentre, single-arm study assessed the safety and performance of the PowerPICC catheters; the studied catheters included PowerPICC, PowerPICC SOLO 2, and PowerGroshong PICC. This study was conducted from 18 June 2020 to 5 January 2022, at 14 centres in 9 European countries (Austria, Belgium, Czech Republic, Denmark, Germany, Italy, Netherlands, Spain, and Switzerland). Good clinical practice principles were adhered to for the design, conduct, recording, and reporting of clinical investigations that

assess the safety and performance of medical devices as described by the International Organization for Standardization (ISO 14155:2011), the principles of the Declaration of Helsinki, applicable sections of national laws and regulations, and EU Medical Device Regulation (Council Regulation 2017/745 of 5 April 2017). This study was registered with the National Institutes of Health clinical trial registry (NCT04263649) (online supplemental table S1). The PowerPICC devices were used for therapy as medically required. Catheter selection for patients, anatomical placement and vessel assessment were performed at the discretion of the clinician. Hence, patients were enrolled non-consecutively. Particular equipment (eg, ultrasound) or assessment protocols (eg, RaPeVa) were used as deemed necessary by the providing clinician during standard of care practice.

Patient population and study devices

Patients were considered for enrolment if they required the use of one of the study devices (PowerPICC, PowerPICC SOLO 2, and PowerGroshong PICC) for either short-term (<30 days) or long-term (≥30 days) vascular access. A key distinction between these devices is their valve-related properties; PowerPICC catheters are open-ended, PowerPICC SOLO 2 catheters have a pressure-activated proximal hub valve, and PowerGroshong catheters are closed-ended (figure 1). Additionally, PowerPICC and PowerPICC SOLO 2 are made of thermoplastic polyurethane, and PowerGroshong PICCs are made of silicone. To be eligible to participate in this study, the PowerPICC catheter selected by the clinician for patient use had to be inserted as standard of care, and the patients were expected to be available for observation through the duration of PICC therapy. Patients were excluded from the study if any of the following criteria applied: (1) the presence of any device-related infection, bacteraemia or septicaemia was known or suspected; (2) insufficient body size to accommodate the size of the implanted device; (3) known or suspected allergic reactions to materials contained in the device; (4) a history of irradiation of the prospective insertion site; (5) previous episodes of venous thrombosis (VT) or vascular surgical procedures occurred at the prospective placement site and (6) presentation of local tissue factors that would have prevented proper device stabilisation and/or

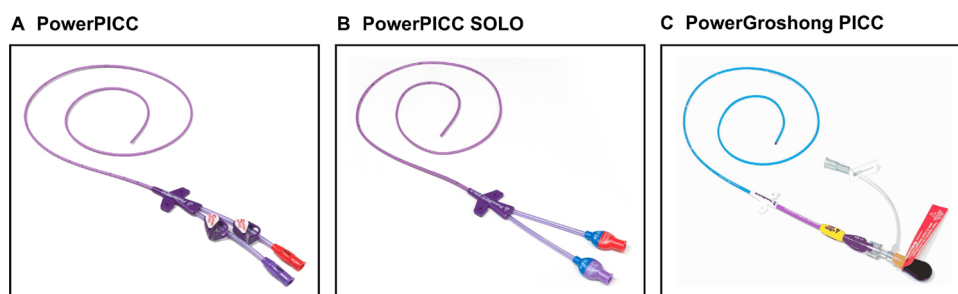


Figure 1 Example images of PowerPICC devices. (A) PowerPICC catheter with open-ended valves. (B) PowerPICC SOLO with pressure-activated proximal hub valve. (C) PowerGroshong PICC with closed-ended, three-position valve which allows fluids to flow in or out but remains closed when not in use. PICC, peripherally inserted central catheter.

access. At Austrian sites only, pregnant and breastfeeding women were excluded from the study, as requested by the Austrian Medical Device Law (MPG). Patients were followed through device removal or 180±14 days postinsertion, whichever came first, in accordance with previous studies.¹³

Clinical study endpoints

Two distinct primary endpoints were assessed in the current study in order to assess both safety and performance outcomes. The primary safety endpoint was the incidence of symptomatic VT with thrombus presence confirmed by ultrasonography or other imaging. The secondary safety endpoints were the incidence of phlebitis, extravasation, vessel laceration (vessel tissue tear), vessel perforation (hole in vessel tissue), local infection, accidental dislodgment, and catheter-related bloodstream infection (CRBSI).

The primary performance endpoint was the proportion of patients who received a PowerPICC device that was successfully used through the completion of therapy. Treatment success was defined as a device that was either still in place at the end of the 180-day observation or had been removed because the therapy was completed (end of useful medical use), changed, or cancelled. Patients who died during the study were considered as having completed the therapy. The secondary performance endpoints were the percent of catheters that were patent (ie, functional) through the completion of therapy, the percent of placement success in a single insertion attempt (based on proper tip location and patency after the first insertion attempt), and usability (based on a post-insertion survey on guidewire/stylet performance). As this was an observational study, no outcome measures were performed outside of what was considered routine practice.

Study design and procedures

The PowerPICC device was used for therapy as medically required, according to the device-specific instructions for use (IFU) and hospital protocols. Information was collected on the study patients and devices before, during, and after catheter insertion while the device was indwelling, and during device removal to assess both safety and performance endpoints. Clinical data for the primary and secondary endpoints were gathered from information collected as part of the patients' routine care, and no study-specific procedures were performed. Choice of the specific PowerPICC type, size, number of lumens, and length required for an individual patient was left to the discretion of the clinician. After the catheter was inserted and procedure data was recorded, the inserting clinician completed a short survey to assess the performance of the guidewire/stylet. The device maintenance (ie, flushing, fluid locking, cleansing) and insertion site maintenance (ie, cleansing, dressing) were performed according to the device-specific IFU, standard medical practice and hospital protocols. The study device was discontinued/

removed when it was deemed medically appropriate by the patient's treating clinician (eg, no longer needed or no longer functioning). Removal of the device was also performed according to the device-specific IFU, standard medical practice, and hospital protocols.

For patients treated in an acute-care setting, the catheter/site was assessed daily until the device was removed or until they were discharged with the device in place (then treated as an outpatient). When the patient was treated in an outpatient setting, the catheter/site was assessed at each therapy visit until the device was removed, or until 180 days postinsertion (±14 days), whichever came first. For the purposes of endpoint calculations, PICCs remaining in place after 180 days were considered to have completed therapy.

Statistical methods

A sample size of 150 participants for each type of device (PowerPICC, PowerPICC SOLO 2 and PowerGroshong PICC) was chosen based on the precision of the point estimates of the primary endpoint (VT), as well as the ability to observe rare adverse events or complications. Assuming the primary endpoint of VT incidence is 10%,¹⁴ the precision of the point estimate is 4.8% with a sample size of 150 (ie, the 95% CI is the point estimate ±4.8%). Additionally, with a sample size of 150, the probability of observing at least one rare complication with a 1% rate is 78%, and the probability of observing at least one rare complication with a 2% rate is 95%.

The incidence rates of symptomatic VT (primary safety endpoint) and of select complications (secondary safety endpoints) and their 95% CI were calculated. The percentage of PICCs remaining in place through the completion of therapy and a 95% CI was calculated based on the total number of PICCs where information about success or failure was documented. The percentage of patent PICCs and the percentage of placement success in a single insertion attempt were calculated together with 95% CIs. The results of the postinsertion user questionnaire were analysed descriptively. No missing value imputation methods were used.

Patient and public involvement

None.

RESULTS

Patient characteristics

A total of 451 patients signed the consent form, of which 1 patient did not meet the eligibility criteria, as illustrated in [figure 2](#). The study device was attempted to be placed in 450 patients (150 patients each for PowerPICC, PowerPICC SOLO 2, and PowerGroshong PICC), and a total of 447 patients received treatment and were included in the evaluable set. PowerPICC SOLO 2 could not be placed in two patients and could not be used for the intended treatment in one patient. Of these patients, 335 were documented as having completed the therapy via

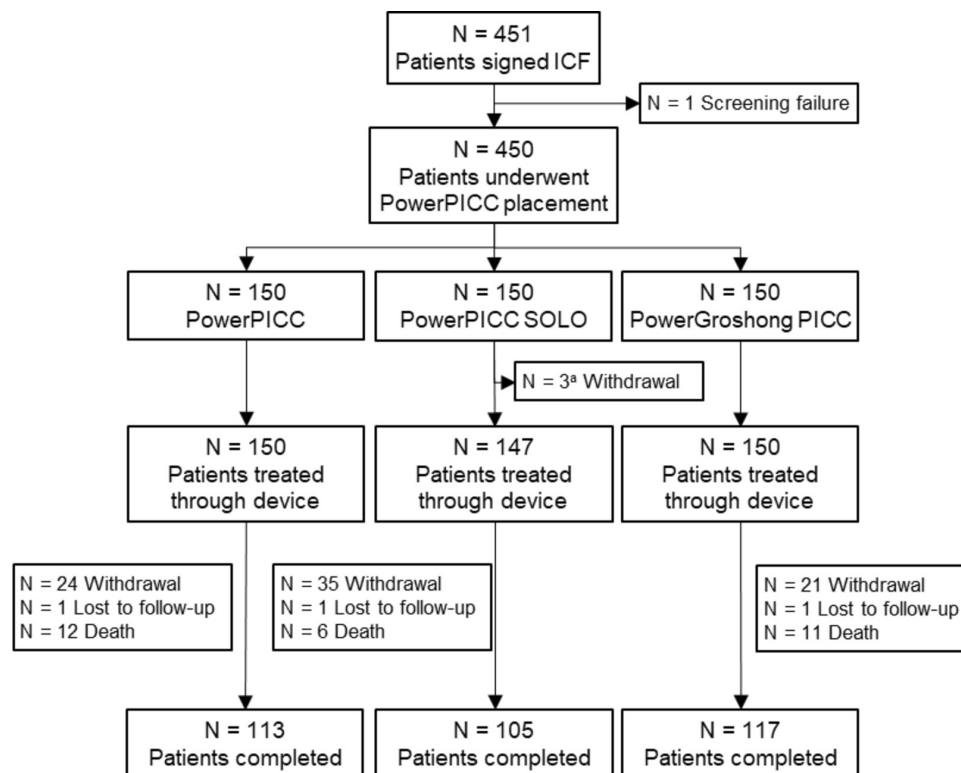


Figure 2 Patient and device disposition of the 29 patients who are documented as being withdrawn because of death, 28 patients had the study device still in place at the time of death and were considered as having completed the therapy for the assessment of the primary performance endpoint and 1 patient was lost to follow-up with no information available about the whereabouts of the device. ^aUnsuccessful placement of study device in two patients; one patient with successful placement of study treatment but could not be treated via the device. ICF, informed consent form; PICC, peripherally inserted central catheter.

the study device (device either removed after therapy or still in place at 180 days), and 115 were documented as non-completers ('withdrawal'). The reasons given for the withdrawal included lost to follow-up in 3 patients, death in 29 patients, adverse events in 53 patients, device deficiency in 19 patients, and 'other' in 11 patients. Catheter to vessel ratio was assessed prior to insertion in 60.4% of all patients.

Of the total patient cohort, 50.2% of patients were women and 49.8% were men. Patient age ranged from 6 to 96 years, with a median of 62.0 years (table 1). The most frequent indication for device placement was the administration of chemotherapy (69.6%). In 78.4% of patients, the PICC was intended for long-term use (ie, for at least 30 days). Device placement was successful in 99.6% of patients.

Safety results

Symptomatic VT after device placement, the primary safety endpoint, was observed in 0.7% (PowerPICC), 1.4% (PowerPICC SOLO), and 2.7% (PowerGroshong PICC) of patients as shown in table 2. When stratified based on the number of lumens, symptomatic VT was observed in six patients (1.6%) among those who received single-lumen PowerPICC catheters and one patient (1.5%) among those who received double-lumen catheters. For secondary safety endpoints, the incidences of the preselected complications across all study devices were 0% for

phlebitis, 0% for extravasation, 4.5% for local infections, 1.6% for CRBSI, and 3.8% for accidental dislodgment. The CRBSI rate per 1000 catheter days was 0.17 across all PowerPICC devices (online supplemental table S2). The device-specific incidence rates for these complications are listed in table 2.

Performance results

For the primary performance outcome, PowerPICC devices remained in place throughout the required therapy time in 81.8% of patients (table 3). For secondary performance outcomes, 93.9% of all PowerPICCs were patent through therapy, and 90.4% could be placed in a single insertion attempt with proper tip location and patency (table 3). When clinicians were asked if the guidewire or stylet facilitated placement of the PICC, the majority of clinicians 'agreed' or 'strongly agreed' that the guidewire and the stylet facilitated placement (94.6% and 94.9%, respectively, online supplemental tables S3 and S4). Use of the guidewire and stylet was reported to be 'easy' or 'very easy' by 93.3% and 94.4% of the clinicians, respectively.

DISCUSSION

This study provides the most robust evidence available on the safety and performance of PowerPICC catheters when used in real-world settings across nine European

Table 1 Patient demographics and baseline characteristics

	PowerPICC N=150	PowerPICC SOLO N=150	PowerGroshong PICC N=150	Any PowerPICC N=450
Sex, N (%)				
Male	80 (53.3)	45 (30.0)	99 (66.0)	224 (49.8)
Female	70 (46.7)	105 (70.0)	51 (34.0)	226 (50.2)
Race, N (%)				
White	141 (94.0)	148 (98.7)	150 (100.0)	439 (97.6)
Asian	4 (2.7)	1 (0.7)	–	5 (1.1)
Black or African American	3 (2.0)	–	–	3 (0.7)
Other	2 (1.3)	1 (0.7)	–	3 (0.7)
Age, median years (min; max)	62.5 (6; 96)	57.5 (21; 92)	63.0 (33; 87)	62.0 (6; 96)
Body height, mean cm (SD)	169.83 (10.61)	169.29 (8.37)	168.62 (8.90)	169.25 (9.34)
Body weight, mean kg (SD)	75.22 (17.51)	77.56 (17.30)	73.45 (13.25)	75.41 (16.20)
BMI, mean (SD)	26.00 (5.54)	27.07 (5.79)	25.83 (4.29)	26.30 (5.27)
Primary diagnosis for current admission, N (%)				
Cardiac	12 (8.0)	–	–	12 (2.7)
Trauma	2 (1.3)	–	–	2 (0.4)
Neurological	7 (4.7)	–	1 (0.7)	8 (1.8)
Infection	33 (22.0)	8 (5.3)	1 (0.7)	42 (9.3)
Surgical	14 (9.3)	14 (9.3)	1 (0.7)	29 (6.4)
Oncological	67 (44.7)	120 (80.0)	141 (94.0)	328 (72.9)
Vascular	4 (2.7)	–	–	4 (0.9)
Gastrointestinal	8 (5.3)	5 (3.3)	5 (3.3)	18 (4.0)
Other	3 (2.0)	3 (2.0)	1 (0.7)	7 (1.6)
Reason for PICC placement, N (%)				
Central venous pressure monitoring	2 (1.3)	–	–	2 (0.4)
Hydration support	12 (8.0)	7 (4.7)	1 (0.7)	20 (4.4)
Intravenous therapy				
Immunoglobulins	1 (0.7)	–	1 (0.7)	2 (0.4)
Antimicrobial/antifungal	56 (37.3)	19 (12.7)	4 (2.7)	79 (17.6)
Blood product and transfusion	13 (8.7)	4 (2.7)	4 (2.7)	21 (4.7)
Chemotherapy	62 (41.3)	116 (77.3)	135 (90.0)	313 (69.6)
Pain management	12 (8.0)	11 (7.3)	7 (4.7)	30 (6.7)
Total parenteral nutrition	22 (14.7)	18 (12.0)	10 (6.7)	50 (11.1)
Vasopressors	4 (2.7)	–	–	4 (0.9)
Limited peripheral access	24 (16.0)	16 (10.7)	85 (56.7)	125 (27.8)
Nutritional support	5 (3.3)	6 (4.0)	4 (2.7)	15 (3.3)
Other	2 (1.3)	–	1 (0.7)	3 (0.7)
Serial blood sampling	29 (19.3)	10 (6.7)	22 (14.7)	61 (13.6)
Serial radiographic studies	1 (0.7)	–	2 (1.3)	3 (0.7)
Thrombocytopenia or coagulopathy	3 (2.0)	–	–	3 (0.7)
Intended access duration, N (%)				
Short term <30 days	67 (44.7)	28 (18.7)	2 (1.3)	97 (21.6)
Long term ≥30 days	83 (55.3)	122 (81.3)	148 (98.7)	353 (78.4)
Number of lumens, N (%)				
Single lumen	94 (62.7)	141 (94.0)	150 (100.0)	385 (85.6)
Double lumen	56 (37.3)	9 (6.0)	–	65 (14.4)

BMI, body mass index; max, maximum; min, minimum; n, number of patients; PICC, peripherally inserted central catheter.

Table 2 Incidence of complications (SAF, N=450, n=447; PowerPICC n=150, PowerPICC SOLO n=147, PowerGroshong n=150)

Complication by device type	Number (%) of patients	95% CI (%)
Primary safety endpoint		
Symptomatic venous thrombosis		
PowerPICC	1 (0.7)	0.0; 3.7
PowerPICC SOLO	2 (1.4)	0.2; 4.8
PowerGroshong	4 (2.7)	0.7; 6.7
Any PowerPICC	7 (1.6)	0.6; 3.2
Secondary safety endpoints		
Phlebitis		
PowerPICC	0 (0.0)	0.0; 2.4
PowerPICC SOLO	0 (0.0)	0.0; 2.5
PowerGroshong	0 (0.0)	0.0; 2.4
Any PowerPICC	0 (0.0)	0.0; 0.8
Extravasation		
PowerPICC	0 (0.0)	0.0; 2.4
PowerPICC SOLO	0 (0.0)	0.0; 2.5
PowerGroshong	0 (0.0)	0.0; 2.4
Any PowerPICC	0 (0.0)	0.0; 0.8
Local infection		
PowerPICC	6 (4.0)	1.5; 8.5
PowerPICC SOLO	13 (8.8)	4.8; 14.6
PowerGroshong	1 (0.7)	0.0; 3.7
Any PowerPICC	20 (4.5)	2.8; 6.8
CRBSI		
PowerPICC	2 (1.3)	0.2; 4.7
PowerPICC SOLO	3 (2.0)	0.4; 5.8
PowerGroshong	2 (1.3)	0.2; 4.7
Any PowerPICC	7 (1.6)	0.6; 3.2
Accidental dislodgment		
PowerPICC	4 (2.7)	0.7; 6.7
PowerPICC SOLO	9 (6.1)	2.8; 11.3
PowerGroshong	4 (2.7)	0.7; 6.7
Any PowerPICC	17 (3.8)	2.2; 6.0
Vein laceration		
PowerPICC	1 (0.7)	0.0; 3.7
PowerPICC SOLO	1 (0.7)	0.0; 3.7
PowerGroshong	0 (0.0)	0.0; 2.4
Any PowerPICC	2 (0.4)	0.1; 1.6
Vein perforation		
PowerPICC	0 (0.0)	0.0; 2.4
PowerPICC SOLO	1 (0.7)	0.0; 3.7
PowerGroshong	0 (0.0)	0.0; 2.4
Any PowerPICC	1 (0.2)	0.0; 1.2

CRBSI, catheter-related bloodstream infection; PICC, peripherally inserted central catheter; SAF, safety analysis set.

countries. Importantly, our prospective, observational, multicentre study approach provides greater generalisability to our findings, as our findings reflect the health

Table 3 Performance and secondary performance endpoints (EVS, N=447, n=445; PowerPICC n=149, PowerPICC SOLO n=146, PowerGroshong PICC n=150)

Primary performance endpoint by catheter type	Number (%) of patients	95% CI (%)
PICCs that remained in place through the required therapy time		
PowerPICC	124 (83.2)	76.2; 88.8
PowerPICC SOLO	111 (76.0)	68.3; 82.7
PowerGroshong	129 (86.0)	79.4; 91.1
Any PowerPICC	364 (81.8)	77.9; 85.3
Secondary performance endpoint	Number (%) of catheters	95% CI (%)
Patent PICCs		
PowerPICC	139 (93.3)	88.0; 96.7
PowerPICC SOLO	135 (92.5)	86.9; 96.2
PowerGroshong	144 (96.0)	91.5; 98.5
Any PowerPICC	418 (93.9)	91.3; 96.0
Placement success in single insertion attempt		
PowerPICC	129 (86.0)	79.4; 91.1
PowerPICC SOLO	129 (86.0)	79.4; 91.1
PowerGroshong	149 (99.3)	96.3; 100.0
Any PowerPICC	407 (90.4)	87.3; 93.0

EVS, evaluable set; PICC, peripherally inserted central catheter.

outcomes that occur when PowerPICC catheters are used in diverse patient populations across European countries with varying standard clinical practices. The patients included in this study ranged from 6 to 96 years of age, and patients were admitted for a variety of health indications, including cardiac, trauma, neurological, oncological, and gastrointestinal conditions. Based on the results of the current study, our findings demonstrate that PowerPICC catheters are safe for treating patients experiencing a variety of health conditions and that these catheters are effectively used by clinicians across Europe.

The incidence of complications was low when using PowerPICC catheters, highlighting the general safety of these devices. For instance, our findings showed that 1.6% of patients developed symptomatic VT across all PowerPICC catheters. This is nearly half that reported by Greene *et al*, which found that of 3790 patients who received a PICC during their hospital stay, ~3.1% of patients developed symptomatic VT.¹⁵ A recent meta-analysis assessing the occurrence of symptomatic deep VT with PICCs revealed an even higher incidence rate of 4.58% for non-intensive care unit (ICU) patients and 5.08% for ICU patients.¹⁶ This observation of higher VT incidence rates with other PICC devices depicts the relatively low incidence of VT that occurs in patients receiving a PowerPICC catheter. In addition, the proportion of patients who developed symptomatic VT was similar between those who received single or double-lumen PowerPICC catheters. Regarding

the study's secondary safety endpoints, there were no observed cases of extravasation or phlebitis and only three cases of vein laceration or vein perforation across all patients who received one of the studied PowerPICC devices. The most common complication observed across all devices was the incidence of local infection, which occurred in 4.5% of all patients. Notably, local infections are also linked to improper catheter placement/maintenance practices (eg, aseptic care) that may be beyond the control of the device.¹⁷ In this regard, our findings are similar to previously published evidence on local infections using PICC devices; Grau *et al* found an incidence of local infections in 4.7% of patients who received PICCs at a single-site hospital in France.¹⁸ Collectively, our findings reveal that complications are rare when using PowerPICC devices, demonstrating that these devices are safely used in the general patient setting in health centres across Europe.

Our study expands on the limited existing knowledge regarding the safety of the PowerPICC catheters. Previous studies were limited to observations of specific patient cohorts, studies that did not distinguish data across specific PowerPICC catheters, or studies that were limited to a single site.^{19–22} Our findings of a 1.6% incidence of symptomatic VT are similar to findings reported by González *et al*, which found an incidence rate of upper extremity deep VT in 2.01% of patients when using PowerPICC catheters. In contrast, the incidence of CRBSI observed in the current study is lower than previously reported values for PowerPICC catheters. Sato *et al* found a CRBSI incidence rate of 1.4 per 1000 catheter days in patients being treated for head and neck cancer when using either PowerGroshong PICC or single-lumen PowerPICC catheters.²¹ Morano *et al* found a CRBSI incidence rate of 0.59 per 1000 catheter days in haematological patients treated at a single haematology centre in Italy when using PowerGroshong or PowerPICC catheters.²³ Our findings of a CRBSI incidence rate of 0.17 per 1000 catheter days indicate that CRBSI occurs less commonly than these previously reported studies focused on specific patient cohorts. Although it is plausible that the differences observed between prior studies and the current findings are reflective of the underlying health conditions of the studied patient populations, our study included patients admitted for a variety of health conditions and reflects the health outcomes that occur in a general hospital setting.

It is important to note that some complications that occur with the use of PICCs may be outside of the control of the device. While the current study's protocols require clinicians to follow care plans according to the IFU of the specific PowerPICC device and hospital standards, medical errors are an unfortunate reality that can negatively impact patient health.²⁴ For instance, bloodstream infections are linked to improper catheter placement/maintenance practices (eg, aseptic care).¹⁷ Successful prevention of complications can also be dependent on the selection of the appropriate catheter (eg, the fewest

number of lumens possible and adequate vein diameter to catheter size ratio based on the patient's needs) and proper catheter management (eg, care of the insertion site, catheter stabilisation, routine flushing).^{16 25 26} Nonetheless, given that there was a low occurrence of complications and that clinicians overwhelmingly reported that PowerPICC devices were easy to use, the current study suggests that these devices are capable of being successfully implemented in the clinical setting and that proper maintenance is readily achievable.

In addition to the safety of these devices, our findings indicate that the PowerPICC catheters perform well in the clinical setting and that clinicians find these devices easy to use. Over 93% of clinicians agreed that both the guidewire and stylet were easy or very easy to use and remove. Indeed, 90.4% of all PowerPICC devices were successfully placed on their first insertion attempt with an overall success rate for placement of 99.6%. The overall success rate across all PowerPICC catheters is higher than compiled values reported by the Society of Interventional Radiology Standard of Practice Committee in their quality improvement guidelines for central venous access, which reports PICC placement success rates of 96% in adult populations.²⁷ Appropriate positioning of PICCs is an important consideration for patient care, as poorly positioned catheter tips may increase the risk of PICC-associated complications including phlebitis, deep vein thrombosis, and catheter malfunctioning.^{28–30} Malpositioned PICCs can also delay the start of patient treatment, require additional hospital resources (eg, repeated PICC procedures and chest X-rays), and increase hospital costs.³¹ Further, receiving multiple PICCs may increase the risk of vascular complications and bloodstream infections.^{32 33} Given that 93.9% of all PowerPICC catheters remained patent through the completion of therapy in the current study, these findings suggest that few patients would require replacement PICCs when using one of the PowerPICC line of catheters. Taken together, these findings highlight the high-performance standard of PowerPICC catheters in the clinical setting, benefitting both clinicians and patients.

Notably, a key distinction between the devices in the PowerPICC catheters is their valve properties; PowerPICC catheters are open-ended, PowerPICC SOLO 2 catheters have a pressure-activated proximal hub valve, and PowerGroshong catheters are closed-ended with a three-position valve (or valves) which allows fluids to flow in or out but remains closed when not in use. Furthermore, the devices differ in their material composition (ie, thermoplastic urethane for PowerPICC and PowerPICC SOLO2; silicone for PowerGroshong). The selection of a PICC with a particular valve type may depend on multiple factors such as patient risk for complications, flow requirements, anticipated PICC maintenance (eg, frequency of flushes), need for intermittent infusion, and considerations of home

discharge. Overall, our findings reveal a similar safety and performance profile across the studied devices indicating that clinicians may select the specific PowerPICC device according to patient needs without concerns of varying efficacy across the PowerPICC devices.

Limitations

There are several limitations of this study. First, standard care practices vary across hospitals and European countries,³⁴ which may have an impact on the safety and performance of PowerPICC catheters. Nonetheless, the low rate of complications and effective performance observed across all PowerPICC devices used at the 14 health centres in the current study indicates that these devices are readily employed in hospital systems across Europe. The specific PowerPICC device used for patients was also not randomised, and PICC selection was performed at the discretion of the clinician in accordance with standard of care practices. Hence, it is plausible that clinician expertise and experience in selecting appropriate VADs may have impacted the outcomes observed in the current study. For instance, the secondary safety outcomes of vein laceration and vein perforation are also inserter-dependent, and only three cases were reported for these outcomes across all observed devices. Inherent differences among catheters may also contribute to findings, such as the availability of catheter sizes; the smallest size available for the PowerGroshong PICC is 5Fr, whereas PowerPICC and PowerPICC SOLO 2 are available in 4Fr sizes. Larger catheter sizes have previously been shown to have an association with an increased incidence of relevant complications (ie, VT).¹⁶ Additionally, although this study included a diverse patient population with patients admitted to the hospital for various health diagnoses, our sample is represented by a large portion of oncology patients. Patients were also admitted for other reasons, including infection, surgical, cardiac, gastrointestinal, and neurological conditions. Nonetheless, the underlying health of the patients included in this study may have impacted the current findings. Due to the relatively low rate of complications, this study did not assess safety outcomes across patients admitted for various diagnoses or different indications for PICC placement. The study design also did not allow for comparisons between other power-injectable compatible devices that may be available, and this approach is a worthy consideration for future studies.

CONCLUSIONS

This study provides robust real-world evidence on the safety and performance of the PowerPICC catheters in diverse patient populations across nine European countries. The low incidence of complications, such as only 1.6% of patients developing symptomatic VT, highlights the safety of these devices. PowerPICC

catheters performed effectively with the majority of catheters successfully placed on their initial attempt and remaining patent through the completion of therapy. These findings support the safe and effective use of PowerPICC catheters in general hospital settings across Europe, providing valuable insights for clinicians and contributing to improved patient care.

Acknowledgements The authors extend their gratitude to the individuals who provided valuable support for data collection: Alfonso Papa, Kim Wedervang, Hans van Overhagen, Ralph Gnannt, Vittorio Fusco, Ulf Teichgräber, Geert Maleux, Rudolf Likar, Fredrik Waneck, Xavier Garcia and Paola Perino. Writing and editorial support were provided by Hayden W. Hyatt and Halit O. Yapici at Boston Strategic Partners (funded by Becton, Dickinson and Company).

Contributors MH, KH and KA designed the study. Data collection was performed by MT, GC, BM, MH, GC, BM, MH, KH and KA contributed to the writing and editing of the manuscript. MH, GC, BM, MH, KH and KA approved of the final version of this manuscript. KA is the guarantor.

Funding This study was funded by Becton, Dickinson and Company as sponsor. All funds have been made payable to all Hospitals involved in the study and not to individual researchers (grant number: N/A).

Competing interests BM consults with Becton, Dickinson and Company (BD), Teleflex, B Braun, IPSA S.r.l. and Medcomp. MH, KH and KA are employees BD and are current equity holders in BD.

Patient and public involvement Patients and/or the public were not involved in the design, or conduct, or reporting, or dissemination plans of this research.

Patient consent for publication Not applicable.

Ethics approval This study involves human participants and was approved by Ethikkommission des Landes Kärnten (MDS-19PICCEU01-AUT-AT01); Ethikkommission der Med. Universität Wien (MDS-19PICCEU01-AUT-AT02); UZ Leuven (MDS-19PICCEU01-BEL-BE01); Kantonale Ethikkommission Zürich (MDS-19PICCEU01-CHE-CH01); Etická komise FN a LF UP Olomouc (MDS-19PICCEU01-CZE-CZ01); Ethikkommission Jena (MDS-19PICCEU01-DEU-DE01); De Videnskabssetiske komiteer Region Syddanmark (MDS-19PICCEU01-DNK-DK01); Comitè d'Ètica d'Investigació Clínica CEIC Girona (MDS-19PICCEU01-ESP-ES01); A.O.R.N. Ospedali dei Colli-Monaldi (MDS-19PICCEU01-ITA-IT01); Comitato etico per la sperimentazione clinica Bolzano (MDS-19PICCEU01-ITA-IT02); Comitato etico interaziendale A.O.U. San Luigi Gonzaga-Torino (MDS-19PICCEU01-ITA-IT03); Comitato Etico Interaziendale A.O.U. Città di Torino (MDS-19PICCEU01-ITA-IT04); Comitato Etico Interaziendale ASO Alessandria (MDS-19PICCEU01-ITA-IT05); Medisch Ethische Toetsingscommissie (METC) (MDS-19PICCEU01-NLD-NL01). Participants gave informed consent to participate in the study before taking part.

Provenance and peer review Not commissioned; externally peer reviewed.

Data availability statement All data relevant to the study are included in the article or uploaded as online supplemental information. No additional data are available.

Supplemental material This content has been supplied by the author(s). It has not been vetted by BMJ Publishing Group Limited (BMJ) and may not have been peer-reviewed. Any opinions or recommendations discussed are solely those of the author(s) and are not endorsed by BMJ. BMJ disclaims all liability and responsibility arising from any reliance placed on the content. Where the content includes any translated material, BMJ does not warrant the accuracy and reliability of the translations (including but not limited to local regulations, clinical guidelines, terminology, drug names and drug dosages), and is not responsible for any error and/or omissions arising from translation and adaptation or otherwise.

Open access This is an open access article distributed in accordance with the Creative Commons Attribution Non Commercial (CC BY-NC 4.0) license, which permits others to distribute, remix, adapt, build upon this work non-commercially, and license their derivative works on different terms, provided the original work is properly cited, appropriate credit is given, any changes made indicated, and the use is non-commercial. See: <http://creativecommons.org/licenses/by-nc/4.0/>.

ORCID iDs

Klaus Hoerauf <http://orcid.org/0000-0001-6600-7190>

Kimberly Alsbrooks <http://orcid.org/0000-0002-8684-7100>

REFERENCES

- 1 Helm RE, Klausner JD, Klemperer JD, *et al.* Accepted but unacceptable: peripheral IV catheter failure. *J Infus Nurs* 2015;38:189–203.
- 2 Carr PJ, Higgins NS, Cooke ML, *et al.* Vascular access specialist teams for device insertion and prevention of failure. *Cochrane Database Syst Rev* 2018;3:CD011429.
- 3 Dychter SS, Gold DA, Carson D, *et al.* Intravenous therapy: a review of complications and economic considerations of peripheral access. *J Infus Nurs* 2012;35:84–91.
- 4 Cheung E, Baerlocher MO, Asch M, *et al.* Venous access: a practical review for 2009. *Can Fam Physician* 2009;55:494–6.
- 5 iData Research. Vascular access devices market size, share & trends analysis | Europe | 2023-2029 | MedSuite | Includes: implantable port market, port needle market, central venous catheter market, and 10 more. 2023. Available: <https://idatasearch.com/product/vascular-access-devices-market-europe/> [Accessed 3 Mar 2024].
- 6 Gorski LA, Hadaway L, Hagle ME, *et al.* Infusion Therapy Standards of Practice, 8th Edition. *J Infus Nurs* 2021;44:S1–224.
- 7 Chopra V, Flanders SA, Saint S, *et al.* The Michigan Appropriateness Guide for Intravenous Catheters (MAGIC): Results From a Multispecialty Panel Using the RAND/UCLA Appropriateness Method. *Ann Intern Med* 2015;163:S1–40.
- 8 Walker G, Todd A. Nurse-led PICC insertion: is it cost effective? *Br J Nurs* 2013;22:S9–15.
- 9 Robinson MK, Mogensen KM, Grudinskas GF, *et al.* Improved care and reduced costs for patients requiring peripherally inserted central catheters: the role of bedside ultrasound and a dedicated team. *JPEN J Parenter Enteral Nutr* 2005;29:374–9.
- 10 Cormode DP, Naha PC, Fayad ZA. Nanoparticle contrast agents for computed tomography: a focus on micelles. *Contrast Media Mol Imaging* 2014;9:37–52.
- 11 Nance JW Jr, Henzler T, Meyer M, *et al.* Optimization of contrast material delivery for dual-energy computed tomography pulmonary angiography in patients with suspected pulmonary embolism. *Invest Radiol* 2012;47:78–84.
- 12 Coyle D, Bloomgarden D, Beres R, *et al.* Power injection of contrast media via peripherally inserted central catheters for CT. *J Vasc Interv Radiol* 2004;15:809–14.
- 13 Campagna S, Gonella S, Berchiolla P, *et al.* A retrospective study of the safety of over 100,000 peripherally-inserted central catheters days for parenteral supportive treatments. *Res Nurs Health* 2019;42:198–204.
- 14 Chopra V, Anand S, Hickner A, *et al.* Risk of venous thromboembolism associated with peripherally inserted central catheters: a systematic review and meta-analysis. *Lancet* 2013;382:311–25.
- 15 Greene MT, Flanders SA, Woller SC, *et al.* The Association Between PICC Use and Venous Thromboembolism in Upper and Lower Extremities. *Am J Med* 2015;128:986–93.
- 16 Bahl A, Alsbrooks K, Gala S, *et al.* Symptomatic Deep Vein Thrombosis Associated With Peripherally Inserted Central Catheters of Different Diameters: A Systematic Review and Meta-Analysis. *Clin Appl Thromb Hemost* 2023;29:10760296221144041.
- 17 Polderman KH, Girbes ARJ. Central venous catheter use. Part 2: infectious complications. *Intensive Care Med* 2002;28:18–28.
- 18 Grau D, Clarivet B, Lotthé A, *et al.* Complications with peripherally inserted central catheters (PICCs) used in hospitalized patients and outpatients: a prospective cohort study. *Antimicrob Resist Infect Control* 2017;6:18.
- 19 Akhtar N, Lee L. Utilization and Complications of Central Venous Access Devices in Oncology Patients. *Curr Oncol* 2021;28:367–77.
- 20 González S, Jiménez P, Saavedra P, *et al.* Five-year outcome of peripherally inserted central catheters in adults: a separated infectious and thrombotic complications analysis. *Infect Control Hosp Epidemiol* 2021;42:833–41.
- 21 Sato R, Kumai T, Hayashi R, *et al.* Peripherally inserted central venous catheters provide safe and easy central venous access in patients with head and neck cancer. *Int J Pract Otolaryngol* 2021;04:e29–33.
- 22 Goossens GA, Grumiaux N, Janssens C, *et al.* SecurAstaP trial: securement with SecurAcath versus StatLock for peripherally inserted central catheters, a randomised open trial. *BMJ Open* 2018;8:e016058.
- 23 Morano SG, Latagliata R, Girmenia C, *et al.* Catheter-associated bloodstream infections and thrombotic risk in hematologic patients with peripherally inserted central catheters (PICC). *Support Care Cancer* 2015;23:3289–95.
- 24 Makary MA, Daniel M. Medical error—the third leading cause of death in the US. *BMJ* 2016;353:i2139.
- 25 Cotogni P, Pittiruti M. Focus on peripherally inserted central catheters in critically ill patients. *World J Crit Care Med* 2014;3:80–94.
- 26 Sharp R, Cummings M, Fielder A, *et al.* The catheter to vein ratio and rates of symptomatic venous thromboembolism in patients with a peripherally inserted central catheter (PICC): a prospective cohort study. *Int J Nurs Stud* 2015;52:677–85.
- 27 Dariushnia SR, Wallace MJ, Siddiqi NH, *et al.* Quality improvement guidelines for central venous access. *J Vasc Interv Radiol* 2010;21:976–81.
- 28 Racadio JM, Doelman DA, Johnson ND, *et al.* Pediatric peripherally inserted central catheters: complication rates related to catheter tip location. *Pediatrics* 2001;107:E28.
- 29 Luciani A, Clement O, Halimi P, *et al.* Catheter-related upper extremity deep venous thrombosis in cancer patients: a prospective study based on Doppler US. *Radiology* 2001;220:655–60.
- 30 Petersen J, Delaney JH, Brakstad MT, *et al.* Silicone venous access devices positioned with their tips high in the superior vena cava are more likely to malfunction. *Am J Surg* 1999;178:38–41.
- 31 Bidgood C. Improving the patient experience with real-time PICC placement confirmation. *Br J Nurs* 2016;25:539–43.
- 32 Yang RY, Moineddin R, Filipescu D, *et al.* Increased complexity and complications associated with multiple peripherally inserted central catheter insertions in children: the tip of the iceberg. *J Vasc Interv Radiol* 2012;23:351–7.
- 33 Kim K, Kim Y, Peck KR. Previous peripherally inserted central catheter (PICC) placement as a risk factor for PICC-associated bloodstream infections. *Am J Infect Control* 2020;48:1166–70.
- 34 Heijink R, Engelfriet P, Rehnberg C, *et al.* A Window on Geographic Variation in Health Care: Insights from EuroHOPE. *Health Econ* 2015;24 Suppl 2:164–77.

Table S1: Ethical Approval numbers for each center

Site number	Institution	Ethics committee that granted study approval	PowerPICC device used
AT-01	Klinikum Klagenfurt	Ethikkommission des Landes Kärnten	PowerPICC SOLO
AT-02	AKH Wien	Ethikkommission der Med.Universität Wien	PowerPICC SOLO
BE-01	University Hospitals KU Leuven	UZ Leuven	PowerPICC SOLO
CH-01	UniversitätsSpital Zürich	Kantonale Ethikkommission Zürich	PowerPICC
CZ-01	Faculty Hospital in Olomouc Czech Republic	Etická komise FN a LF UP Olomouc	PowerPICC PowerPICC SOLO
DE-01	Universitätsklinikum Jena	Ethikkommission Jena	PowerPICC
DK-01	Sygehus Sønderjylland, Sønderborg Denmark	De Videnskabetiske komiteer Region Syddanmark	PowerPICC SOLO
ES-01	Hospital Universitari de Girona	Comitè d'Ètica d'Investigació Clínica CEIC Girona	PowerPICC
IT-01	Azienda Ospedaliera dei Colli-Monaldi Napoli Italy	A.O.R.N. Ospedali dei Colli - Monaldi	PowerGroshong
IT-02	Ospedale Centrale di Bolzano	Comitato etico per la sperimentazione clinica Bolzano	PowerPICC PowerPICC SOLO
IT-03	Civil Hospital Ivrea	Comitato etico interaziendale A.O.U. San Luigi Gonzaga - Torino	PowerGroshong
IT-04	University Hospital Molinette, Torino	Comitato Etico Interaziendale A.O.U. Città di Torino	PowerPICC PowerGroshong

IT-05	Azienda Ospedaliera Nazionale SS. Antonio e Biagio e Cesare Arrigo	Comitato Etico Interaziendale ASO Alessandria	PowerGroshong
NL-01	HagaZiekenhuis, Den Haag	Medisch Ethische Toetsingscommissie (METC)	PowerPICC

Table S2: Number of catheter related bloodstream infections per 1000 PICC days (SAF, N = 450)

Type of device	Number of patients with catheter related bloodstream infection	Sum of device lifetimes [days]	Number of bloodstream infections per 1000 PICC days
PowerPICC (N = 150)	2	6746.6	0.30
PowerPICC SOLO (N = 150)	3	12989	0.23
PowerGroshong PICC (N = 150)	2	21054.2	0.09
Any Power PICC (N = 450)	7	40789.8	0.17

N = number of patients in data set, PICC = Peripherally Inserted Central Catheter, SAF = safety analysis set.

Table S3: Usability: ease of insertion - guidewire (SAF, N = 450)

	Number (%) of patients			
	PowerPICC N = 150	PowerPICC SOLO N = 150	PowerGroshong PICC N = 150	Any Power PICC N = 450
Difficulty to use the guidewire				
Very difficult	-	-	-	-
Difficult	2 (1.3)	1 (0.7)	1 (0.7)	4 (0.9)
Neither difficult nor easy	5 (3.3)	6 (4.0)	15 (10.0)	26 (5.8)
Easy	64 (42.7)	25 (16.8)	77 (51.3)	166 (37.0)
Very easy	79 (52.7)	117 (78.5)	57 (38.0)	253 (56.3)
Force to advance the guidewire acceptable				
No	-	-	-	-
Yes	149 (100.0)	149 (100.0)	150 (100.0)	448 (100.0)
Difficulty to withdraw/remove the guidewire				
Very difficult	-	-	-	-
Difficult	-	-	-	-
Neither difficult nor easy	3 (2.0)	2 (1.3)	18 (12.0)	23 (5.1)
Easy	64 (43.0)	24 (16.1)	73 (48.7)	161 (35.9)
Very easy	82 (55.0)	123 (82.6)	59 (39.3)	264 (58.9)
Force to withdraw/remove the guidewire acceptable				
No	1 (0.7)	-	-	1 (0.2)
Yes	149 (99.3)	149 (100.0)	150 (100.0)	448 (99.8)
Agreement on integrity				
No	1 (0.7)	-	-	1 (0.2)
Yes	149 (99.3)	149 (100.0)	150 (100.0)	448 (99.8)
Guidewire facilitated placement				
Strongly disagree	-	-	-	-
Disagree	1 (0.7)	-	-	1 (0.2)
Neither agree nor disagree	3 (2.0)	4 (2.7)	16 (10.7)	23 (5.1)
Agree	79 (52.7)	26 (17.4)	72 (48.0)	177 (39.4)
Strongly agree	67 (44.7)	119 (79.9)	62 (41.3)	248 (55.2)

Note: N = 0 is shown as '-'.
 N = number of patients, PICC = Peripherally Inserted Central Catheter, SAF = safety analysis set.

Table S4: Usability: ease of insertion - stylet (SAF, N = 450)

	Number (%) of patients			
	PowerPICC N = 150	PowerPICC SOLO N = 150	PowerGroshong PICC N = 150	Any Power PICC N = 450
Difficulty to use the stylet				
Very difficult	-	-	-	-
Difficult	-	-	-	-
Neither difficult nor easy	1 (1.6)	3 (2.1)	16 (10.7)	20 (5.6)
Easy	48 (77.4)	27 (18.8)	78 (52.3)	153 (43.1)
Very easy	13 (21.0)	114 (79.2)	55 (36.9)	182 (51.3)
Force to advance the stylet/catheter acceptable				
No	-	1 (0.7)	-	1 (0.3)
Yes	61 (100.0)	141 (99.3)	150 (100.0)	352 (99.7)
Difficulty to withdraw/remove the stylet/catheter				
Very difficult	-	-	-	-
Difficult	1 (1.7)	-	-	1 (0.3)
Neither difficult nor easy	-	1 (0.7)	8 (5.4)	9 (2.5)
Easy	46 (76.7)	25 (17.4)	83 (55.7)	154 (43.6)
Very easy	13 (21.7)	118 (81.9)	58 (38.9)	189 (53.5)
Force to withdraw/remove the stylet/catheter acceptable				
No	1 (1.6)	-	-	1 (0.3)
Yes	61 (98.4)	143 (100.0)	149 (100.0)	353 (99.7)
Agreement on integrity				
No	-	-	-	-
Yes	62 (100.0)	144 (100.0)	150 (100.0)	356 (100.0)
Stylet/catheter facilitated placement				
Strongly disagree	-	-	-	-
Disagree	-	-	-	-
Neither agree nor disagree	1 (1.6)	3 (2.1)	14 (9.3)	18 (5.1)
Agree	44 (71.0)	29 (20.1)	78 (52.0)	151 (42.4)
Strongly agree	17 (27.4)	112 (77.8)	58 (38.7)	187 (52.5)
Sherlock 3CG pre-loaded stylet question^a answered				
Yes	48 (32.0)	85 (57.0)	125 (83.3)	258 (57.5)
Not applicable/not used	102 (68.0)	64 (43.0)	25 (16.7)	191 (42.5)
Sherlock 3CG pre-loaded stylet aided				
Strongly disagree	-	-	-	-
Disagree	-	-	-	-
Neither agree nor disagree	-	-	14 (11.2)	14 (5.4)
Agree	31 (64.6)	25 (29.4)	75 (60.0)	131 (50.8)

Strongly agree	17 (35.4)	60 (70.6)	36 (28.8)	113 (43.8)
----------------	-----------	-----------	-----------	------------

Note: N = 0 is shown as '-'.^a

^a Was a Sherlock 3CG used for placement?

N = number of patients, PICC = Peripherally Inserted Central Catheter, SAF = safety analysis set.