

AperTO - Archivio Istituzionale Open Access dell'Università di Torino

Comprehensive two-dimensional gas chromatography and food sensory properties: potential and challenges

This is the author's manuscript

Original Citation:

Availability:

This version is available <http://hdl.handle.net/2318/158481> since 2016-12-01T13:07:17Z

Published version:

DOI:10.1007/s00216-014-8248-z

Terms of use:

Open Access

Anyone can freely access the full text of works made available as "Open Access". Works made available under a Creative Commons license can be used according to the terms and conditions of said license. Use of all other works requires consent of the right holder (author or publisher) if not exempted from copyright protection by the applicable law.

(Article begins on next page)



UNIVERSITÀ DEGLI STUDI DI TORINO

The final publication is available at Springer via <http://dx.doi.org/10.1007/s00216-014-8248-z>

Two-dimensional comprehensive gas chromatography and food sensory properties: potentials and challenges

Chiara Cordero^{1*}, Johannes Kiefl², Peter Schieberle³, Stephen E. Reichenbach⁴ and Carlo Bicchi¹

Authors' affiliation:

¹*Dipartimento di Scienza e Tecnologia del Farmaco, Università degli Studi di Torino,
Via Pietro Giuria 9, I-10125 Torino, Italy*

²*Research & Technology, Symrise AG
Muehlenfeldstrasse 1, 37603 Holzminden, Germany*

³*Chair for Food Chemistry, Technical University of Munich;
Lise-Meitner-Strasse 34, 85354 Freising, Germany*

⁴*Computer Science and Engineering Department, University of Nebraska – Lincoln
260 Avery Hall, Lincoln, NE 68588-0115, USA*

* Address for correspondence:

*Dr. Chiara Cordero - Dipartimento di Scienza e Tecnologia del Farmaco, Università di Torino, Via Pietro
Giuria 9, I-10125 Torino, Italy – e-mail: chiara.cordero@unito.it ; phone: +39 011 6707662; fax: +39 011
2367662*

25 **Abstract**

26 Modern -omics disciplines dealing with food flavor focus the analytical efforts on the elucidation of
27 sensory-active compounds including all possible stimuli of multimodal perception (aroma, taste, texture,
28 etc.) by means of a comprehensive, integrated treatment of sample constituents, such as physicochemical
29 properties, concentration in-the-matrix, and sensory properties (odor/taste quality, perception threshold).
30 Such analyses require detailed profiling of known bioactive components as well as advanced fingerprinting
31 techniques to catalog sample constituents comprehensively, quantitatively, and comparably across samples
32 Multidimensional analytical platforms support comprehensive investigations required for flavor analysis by
33 combining information on analytes' identities, physicochemical behaviors (volatility, polarity, partition
34 coefficient, and solubility), concentration, and odor quality. Interestingly, unlike other *omics*, flavor
35 metabolomics and sensomics include the final output of the biological phenomenon (i.e. sensory
36 perceptions) as an additional analytical dimension, which is specifically and exclusively triggered by the
37 analyzed chemicals. However, advanced -omics platforms, which are multidimensional by definition, pose
38 challenging issues not only in terms of hyphenation with detection systems and sample preparation, but
39 also in data elaboration and processing. The large number of variables collected during each analytical run
40 provides a high level of information, but requires appropriate strategies to fully exploit this potential.
41 This review focuses on advances in two-dimensional comprehensive gas chromatography (GC×GC) and
42 analytical platforms combining GC×GC with olfactometry, chemometrics, and quantitative assays for food
43 sensory analysis to assess the quality of a given product. Sections review instrumental advances and
44 hyphenations, automation in sample preparation, data elaboration, and a selection of applications.

45
46
47

48 **Key-words:**

49 Two-dimensional Comprehensive Gas Chromatography; Gas Chromatography - Olfactometry; Sensomics;
50 Food aroma; Data Processing; Non-targeted analysis; Targeted profiling; Sample preparation; Solid Phase
51 Microextraction; Headspace Analysis; High Concentration Capacity Headspace techniques; Fingerprinting;
52 Multidimensional Gas Chromatography

53

54 **Introduction**

55 *Targeted omics for food sensory quality objectification*

56 Modern *-omics* disciplines dealing with food quality or authentication (foodomics, flavour metabolomics,
57 flavoromics, sensomics [1-5]) investigate sample constituents considering collectively primary and
58 secondary metabolites, and compounds generated or modified by e.g., thermal treatments and/or
59 enzymatic activity, processing, storage, and/or biotechnological treatments.

60 Sensomics and flavoromics, in particular, focus the analytical efforts on elucidating sensory-active
61 compounds and on all possible stimuli of multimodal perception (aroma, taste, texture, etc.) by means of a
62 comprehensive, integrated treatment of sample constituents and their related attributes, such as
63 physicochemical properties, concentration in-the-matrix, or sensory properties [4]. Such analyses require
64 detailed profiling of known bioactive components as well as advanced fingerprinting techniques to catalog
65 sample constituents comprehensively, quantitatively, and comparably across samples [1,6].

66 Conventional, well-established approaches adopted in *-omics* studies for food aroma characterization [7]
67 aim to isolate, identify and quantify key-aroma compounds by combining extraction (liquid-liquid extraction
68 or more effective processes such as Solvent Assisted Flavor Evaporation (SAFE), Simultaneous Distillation
69 Extraction (SDE), Solid Phase Extraction (SPE), and Supercritical fluids Extraction (SFE)) odorant detection by
70 GC-Olfactometry (GC-O), identification, and subsequent accurate quantitation. These approaches are not
71 only fundamental to describe flavor composition and key-components but also for high-throughput
72 screenings and fingerprinting [8].

73 This review focuses on advances in two-dimensional comprehensive gas chromatography (GC×GC) and
74 analytical platforms combining GC×GC with olfactometry, chemometrics, and quantitative assays for food
75 sensory quality assessment. Sections review instrumental advances and hyphenations, automation in
76 sample preparation, and a selection of applications. A section also is devoted to bi-dimensional data
77 elaboration, being that this step of the analytical process is fundamental to fully exploiting all the
78 information included in each analytical run.

79 *The key-role of multidimensionality in food aroma investigations*

80 Food aroma perception is a complex biological phenomenon triggered by certain volatile molecules, mostly
81 hydrophobic, sometimes occurring in trace-level concentrations (at mg/Kg or µg/Kg levels). These
82 molecules must be able to interact with a complex array of Odorant Receptors (ORs) expressed by Olfactory
83 Sensory Neurons (OSNs) in the olfactory epithelium [9-12]. Perception is triggered by specific ligand-
84 receptor interactions and the simultaneous activation of different ORs generates a complex pattern of
85 signals (i.e., the Receptor Code) that is subsequently integrated by the peripheral and central nervous
86 system. Thus, an accurate and comprehensive chemical characterization of the mixture of potential ligands
87 (i.e., the Chemical Odor Code) is fundamental to (i) understand what drives olfactory perception and (ii) to
88 objectify food aroma evaluation.

89 From this perspective, analytical chemistry and separation science play important roles in basic studies of
90 flavor chemistry [13], and modern multidimensional analytical (MDA) platforms are valuable tools for this
91 intriguing field [14-16].

92 Multidimensional platforms support the comprehensive investigations required for flavor chemistry
93 research by combining information on: (a) analytes' identities provided by Mass Spectrometry (MS)
94 through exact mass assignment (High-Resolution Mass Spectrometry - HRMS), diagnostic fragmentation
95 patterns provided by Electron Impact ionization (EI-MS), and/or multiple reaction monitoring by tandem
96 MS techniques (MS/MS or MSⁿ); (b) analytes' physicochemical behaviors based on volatility, polarity,
97 partition coefficient, and solubility; (c) analytes' quantitation based on true (absolute) concentrations in
98 samples or relative abundances, and last, but not least, (d) analytes' odor quality. The latter is possible by
99 implementing olfactometric detection, i.e., human assessors detect odor-active compounds as they elute
100 from a GC column [17-19].

101 Flavor metabolomics and sensomics, unlike other *omics*, include, as an additional analytical dimension, the
102 final output of the biological phenomenon, i.e., sensory perception, which is specifically and exclusively
103 triggered by the analyzed chemicals. However, advanced *-omics* platforms, which are multidimensional by
104 definition, pose challenging issues not only in terms of hyphenation with detection systems and sample
105 preparation, but also in data elaboration and processing. The large number of variables collected during
106 each analytical run provides a high level of information, but requires appropriate strategies to fully exploit
107 this potential.

108

109 *Emerging fields for the application of GC×GC*

110 Apart from quality control aspects, the focus in food aroma analysis over the last five decades has moved
111 from characterizing key odorants and their formation in food to understanding the interaction-relationship
112 with flavor perception, personal behavior, and health. Although most of the key odorants of commonly
113 known foods have been identified [12], more complex questions remain; for example, the role of odorants
114 in odor and flavor perception is poorly understood. One way to study such interactions is to correlate the
115 Chemical Odor Code with sensory data and extract those relevant odorants that modulate the different
116 flavor sensations of a given food [20]. However, for this purpose, ideally, the entire set of key odorants
117 should be measured without discriminating between the highly abundant and chromatographically well-
118 resolved peaks. GC×GC has proven to be a valuable tool to quickly perform a comprehensive assessment of
119 such odorants.

120 Dunkel et al. [12] showed that the development of the Chemical Odor Code of foods is strongly influenced
121 by the manufacturing process, giving highly connected key odorant patterns. Although the development of
122 analytical methods based on GC×GC can be more time consuming than with 1D GC [21], GC×GC based

platforms have greater capacity to resolve such key odorant patterns from different foods, leading to a more effective profiling.

The odorants in common, so-called “generalists”, frequently occur in fermented, aqueous thermally processed (boiled, cooked), and thermally processed (roasted, deep-fried, baked) foods, and are generated from carbohydrates, amino acids, and unsaturated fatty acids as ubiquitous biosynthetic precursors [12]. However, many key odorants are “individualists”, which are unique to certain foods, and so analysis with 1D GC requires optimization of individual methods for proper analytical characterization, if several foods are analyzed routinely. Thus, gas chromatographic analysis of key odorants benefits from enhanced peak capacity of GC×GC and mass spectrometric capabilities to assess large sample sets for: (a) statistical correlations, e.g., with sensory data, and (b) faster characterization of the Chemical Odor Code in different foods.

Advances in two-dimensional comprehensive gas chromatographic analytical platforms

Sensomics aims [4, p.417] “to map the combinatorial code of aroma and taste-active key molecules, which are sensed by human chemosensory receptors and are then integrated by the brain...”. Methods for aroma characterization involve several key steps: (a) extraction and isolation of volatiles; (b) concentration of extracts; (c) pre-separation and fractionation to reduce sample dimensionality [22]; (d) chromatographic separation and selection of intense aroma compounds; (e) identification of odor active compounds and other sample/fraction major components; (f) quantitation; and (g) validation of the aroma contributions by recombination and omission experiments [4,7,12,13].

Some of these discrete and time-consuming steps can be merged and combined in one single analytical system, e.g., GC×GC platforms that take advantage of the vast experiences and instrumental solutions already available for Multidimensional Gas Chromatography (MDGC) [14-16,19,23].

MDGC plays an important role in flavor research, which often requires in-depth investigations [14,15,24], and has a long history, although its widespread application still remains unfulfilled after many years [15,22]. The driving force behind the development of MDGC in the early days of capillary GC was the recognition that, for complex samples, single-columns were often inadequate to provide the expected analytical results. The demand for resolved chromatographic peaks was the force behind this search that resulted in the first instrumental arrangement for comprehensive two-dimensional GC separations in the early 1990s [26,27].

A single GC column has a theoretical informing power (or peak capacity) of about 500-600 (i.e., 500-600 evenly distributed peaks (compounds) can be separated in a single analysis [28]); however, peaks are neither evenly nor randomly distributed in a chromatogram because of sample dimensionality, i.e., the degree of chemical correlation among analytes/constituents. This is particularly true for food samples of vegetable origin whose volatile fraction is characterized by secondary metabolites with common/similar

moieties because of their common biosynthetic pathways. On the other hand, the complex pattern of volatiles produced by thermal processing of food (e.g., roasted coffee or nuts) also creates separation challenges due to the high number of structurally correlated analytes formed from common precursors, for example, homologues and isomers of alkenes, aldehydes, ketones, alcohols, acids, esters, lactones, and phenols, and series of heterocyclic compounds such as furans, pyrazines pyrroles, thiophenes, pyridines, thiazoles, and oxazoles. As a consequence, the required system peak capacity must be much higher than the actual number of compounds in the sample to achieve complete resolution. The result of these factors is that complex samples have a high likelihood of multiple peak co-elutions in a single separation and, according to Davis and Giddings [22,29,30], may require multidimensional separations.

With this perspective, it was immediately evident that GC×GC provides substantial advantages for the detailed characterization of complex mixtures such as some food-derived volatiles, including odor-active compounds responsible for sensory attributes.

One of the first applications in this field was presented by Adahchour et al. [31,32], who investigated the informative potential of GC×GC-TOFMS for the detailed analysis of extracts from milk-derived products (dairy and non-dairy sour cream and dairy spread). The analytical platform was equipped with a Longitudinally Modulated Cryogenic System (LMCS) and consisted of a ^1D 15 m × 0.25 mm ID × 0.25 μm d_f CP-Sil 5 CB low bleed/MS phase column (Varian-Chrompack, Middelburg, the Netherlands) connected, via press-fit, to a ^2D 0.8 m × 0.1 mm ID × 0.1 μm d_f BPX-50 column (SGE Europe, Milton Keynes, UK). Extracts, obtained with well-established techniques, i.e., Solvent Assisted Flavor Evaporation (SAFE) and Cold Finger (CF) distillation [33,34], were analyzed under optimized separation conditions to fully exploit the system's potential. As stated by the authors, the results convincingly showed the merits of the technique for both the overall qualitative characterization (detailed profiling) of volatiles from milk-derived products and the targeted key-flavor components quantitation. Compared to 1D-GC-TOFMS, the quality of the mass spectra obtained after GC×GC separation was higher and enabled more reliable identifications, especially for those analytes closely eluting with interfering matrix compounds. The enhanced overall chromatographic resolution also facilitated quantitation of target compounds, such as methional and sotolon, that were found to be present in the extracts at mg/kg concentration, whereas 1D-GC-TOFMS gave a 100-fold overestimation. The need for further improvements of the technique by devising alternative separation strategies, as reported by authors in the concluding remarks, were the seeds of the subsequent instrumental developments that appeared a few years later.

This discussion of the advances of the analytical platform would be incomplete without a brief discussion of the GC×GC core component, i.e., the modulator. The characterization of key odorants requires effective trapping and release of highly volatile analytes, most of them being responsible for distinct odor notes of some food products and some present in trace amounts. To obtain a suitable band-focusing before entering in the ^2D column, while avoiding breakthrough, dual-stage thermal modulators with cooling media

193 (CO₂ or liquid N₂) have been prevalent. They also enable narrow-bore ²D columns that improve the Signal-
194 to-Noise ratio (SNR) [35] and thus the overall sensitivity of a method. Only a few recent studies have been
195 conducted with flow modulators (FM) and/or cryogenic-free thermal modulators, but, in the authors'
196 opinion, they are worthy of note because they may facilitate adoption of this technique in food quality-
197 control laboratories.

198 Manzano et al. [36] recently studied the volatile fraction of roasted almonds using a commercial flow-
199 modulator by Agilent (Little Falls, DE, US), based on the Capillary Flow Technology (CFT). The authors
200 applied Static Headspace extraction (S-HS) on raw and roasted almonds (*Prunus amygdalus* L. var. *dulcis*) of
201 the Spanish cultivar *Largueta*, and tested different column stationary phase combinations to obtain
202 informative separation patterns. The system was equipped with a Flame Ionization Detector (FID) and
203 analyte identification with references was limited to 43 targets. Although this study is interesting with
204 respect to the potential adoption of a simpler and cost-effective modulator, its main limitation is the
205 absence of mass spectrometric detection, thus limiting the investigation to external standard confirmable
206 analytes and/or to fingerprinting classifications.

207 A study by Tranchida et al. [37], presenting a flexible loop-type flow modulator for GC×GC-FID, discussed its
208 potential for the detailed characterization of spearmint essential oil. The interface consisted of a self-made
209 capillary flow modulator with seven ports connected to an auxiliary pressure source via two branches, to
210 the ¹D and ²D, to a waste branch and a variable modulation loop (2 ports). The spearmint essential oil was
211 separated on a ¹D Enantio Selective (ES) stationary phase coated column, a MEGA-DEX DET-Beta (2,3-
212 diethyl-6-tert-butyl dimethylsilyl-β-CD) 20 m × 0.10 mm ID × 0.10 μm *d_f* (MEGA, Legnano, Italy), coupled to
213 a ²D Supelcowax-10 (polyethylene glycol) 2.5 m × 0.25 mm ID × 0.25 μm *d_f* (Supelco, Bellefonte, PA, USA).
214 Although a satisfactory separation was achieved, the authors stated that further research was necessary to:
215 (a) improve the transfer system to generate well shaped peaks and (b) obtain close-to-optimum second-
216 dimension velocities while keeping an adequate overall sensitivity. More recently, the same research group
217 [38] presented improvements of the flexible loop-type FM, with which citrus essential oil components were
218 effectively separated without a remarkable loss of sensitivity by varying capillary-loop capacity. In this
219 study, tandem MS detection with a triple quadrupole system was employed.

220 Instrumental advances on GC×GC platforms that implement most of the well established techniques of the
221 flavor chemistry community have been defined by Marriott and co-workers as “multi-multidimensional”
222 approaches [39]. In 2010, Maikhunthod et al. [40] presented an instrumental solution that allowed
223 switching between comprehensive two-dimensional gas chromatography and targeted multidimensional
224 gas chromatography system (i.e., switchable GC×GC/targeted MDGC). A schematic diagram of the system is
225 shown in **Figure 1**. The system enabled separate and independent analyses by 1D-GC, GC×GC, and targeted
226 MDGC with the additional possibility of switching from GC×GC to targeted MDGC any number of times
227 throughout a single analysis. Using a Deans switch microfluidics transfer module and a cryotrapping device

(CT), the ¹D column effluent could be directed to either one of two ²D columns in a classical heart-cutting operation. The function of the CT was to focus effectively and rapidly remobilize solute bands to the respective second columns. A short second column enabled GC×GC operation and a longer column was used for targeted MDGC. The system's operational performance parameters were validated by using a mixture of volatiles of interest in the flavor and fragrance field, and on lavender essential oil. Figures of merit were related mainly to obtaining better resolved peaks by a targeted separation on a longer ²D column by diverging specific regions of a GC×GC separation in which co-elutions occurred. Co-elutions in fact prevent reliable identification and quantitation of target analytes.

The potentials of hyphenated and multi-multidimensional systems to study aroma-impact compounds were exploited by Chin et al. [41] in a study focused on coffee brews and Australian wines (Merlot and a blend of Sauvignon Blanc and Semillon). The authors implemented a system capable of GC-Olfactometry (GC-O) and GC×GC with various detectors (TOFMS, FID, and FPD in sulfur mode). In aroma screening modality, the system utilized a ¹D column (DB-FFAP; 15 m × 0.25 mm ID × 0.25 μm *d_f*) connected by means of a Y-split union press fit to a deactivated fused silica tubing (55 cm × 0.1 mm ID) to transfer one half of the effluent to the olfactory port. The other outlet directed the remaining flow to a ²D column (DB-5; 1.1 m × 0.1 mm ID × 0.1 μm *d_f*) connected to a FID detector. A thermal modulator (LMCS) was installed after the Y-split union along the head of the ²D column. The detection frequency method Surface of Nasal Impact Frequency (GC-SNIF) [42] was used for GC-O screening of the volatiles isolated by Solid Phase Extraction (SPE).

Several character-impact odorants were tentatively identified by correlating data of GC×GC-FPD with TOFMS. In particular, the most odor-active analytes from coffee SPE extracts were reported to be 2-methyl-2-butenal, 2-(methoxymethyl)-furan, dimethyl trisulfide, 2-ethyl-5-methyl-pyrazine, 2-octenal, 2-furancarboxaldehyde, 3-mercapto-3-methyl-1-butanol, 2-methoxy-3-(2-methylpropyl)-pyrazine, 2-furanmethanol, and isovaleric acid. From the Australian wines, some varietals' aroma compounds were also identified: 1-octen-3-ol, butanoic acid, and 2-methylbutanoic acid were present in both the Merlot and the Sauvignon Blanc + Semillon (SV) blend with high aroma potency. On the other hand, several co-eluting peaks of ethyl 4-oxo-pentanoate, 3,7-dimethyl-1,5,7-octatrien-3-ol, (Z)-2-octen-1-ol, and 5-hydroxy-2-methyl-1,3-dioxane were suggested to contribute to the Merlot wine aroma; whereas (Z)-3-hexen-1-ol, β-phenylethyl acetate, hexanoic acid, and co-eluting peaks of 3-ethoxy-1-propanol and hexyl formate contributed to the SV wine aroma character. Of the volatile sulfur compounds, 2-mercapto-ethyl acetate was found to add a *fruity, brothy, meaty, sulfur* odor to the Australian wine aroma. The approach of integrating GC-O with concurrent GC×GC analysis successfully revealed the wide range of volatiles present within the most informative odor regions of the 2D chromatograms. The correlation across various GC×GC modalities, coupled to MS identification and sulfur specific detection, provided selective and compound-specific detection to support identification.

262 A further advancement of this platform was presented recently by the same authors [40]. The newer
263 system was capable of performing 1D-GC, GC×GC, and targeted heart-cut MDGC (H/C MDGC) using
264 olfactometry (O), FID, and/or quadrupole MS detection. The system was equipped with a liquid carbon
265 dioxide cryotrapping device (CT) for multiple SPME desorption [43] and H/C MDGC, an olfactory port, a
266 Deans switch (DS) device, a 2-way effluent splitter (ES) based on microfluidics technology, and a thermal
267 modulator (Everest model Longitudinally Modulated Cryogenic System - LMCS). The final configuration is
268 shown in **Figure 2**.

269 The column configuration was: ¹D column DB-FFAP (30 m × 0.25 mm ID × 0.25 μm *d_f*), ²D GC×GC short
270 column for BPX5 (0.9 m × 0.10 mm ID × 0.10 μm *d_f*), and ²D MDGC long column DB-5 MS column (30 m ×
271 0.25 mm ID × 0.25 μm *d_f*). The effluent from the short ²D column outlet was split equally to a FID and the
272 olfactory port by a Y-type device and two deactivated fused silica capillaries (55 cm × 0.10 mm ID). The
273 effluent from the long ²D column outlet was split by the ES device in a ratio of 1:1 and directed to the MS
274 detector via a transfer line (80 cm × 0.10 mm ID) heated at 240 °C and the olfactory port via another
275 transfer line (75 cm × 0.10 mm ID).

276 The integrated analytical system enabled an investigation strategy combining GC×GC-FID/GC-O for an initial
277 screening of odor regions to identify target odor regions (GC-O) and a rapid qualitative and quantitative
278 profiling of the entire complex mixture (GC×GC-FID). The H/C MDGC provided a better separation of
279 targeted regions, depending on the combination of the selected stationary phases, and a contemporary
280 qualification of odor quality/intensity accompanied by analyte identification by qMS (MDGC-MS/O).

281 Experimental results on Shiraz wine volatiles demonstrated the effectiveness of the hyphenated platform,
282 enabling the tentative identification of some odorants: acetic acid, octen-3-ol, and ethyl octanoate as
283 relevant aroma contributors and the determination of β-damascenone (*floral* odor) well separated from
284 hexanoic acid (*sweaty* odor). An analysis of dried spices [44] also indicated the usefulness of the approach
285 by successful identification of character impact-odorant changes during shelf-life. With the integrated
286 system for GC×GC-FID/GC-O combined with automated Headspace Solid Phase Microextraction (HS-SPME),
287 some monoterpenoids were positively correlated with the freshness of the fennel samples; with β-pinene,
288 sabinene, β-myrcene, α-phellandrene, and neo-allo-ocimene found to be more abundant in fresh samples
289 than in five year old products.

290 Recently, Mommers et al. [45] proposed a tunable secondary dimension selectivity system for GC×GC-MS.
291 The tunable system consisted of three capillary columns, different in terms of selectivity and retention
292 mechanisms, one installed as the primary column (¹D) and two, serially coupled, as the secondary column
293 (²D). The ¹D column was a 30 m × 0.25 mm ID, 1 μm *d_f* VF1 MS column (100% dimethyl polysiloxane) and
294 the ²D second dimension consisted of two columns coupled in series: a polar 1 m × 0.1 mm ID, 0.1 μm *d_f*
295 Wax-HT® (100% polyethylene glycol) and a medium polarity 2 m × 0.1 mm ID, 0.2 μm *d_f* VF17 MS column
296 (50% phenyl 50% dimethylpolysiloxane). The contribution of the first of the ²D columns was varied by

altering its effective length, by sliding it stepwise back or forward through the modulator and/or by applying a temperature offset with respect to the main oven. By adjusting the contribution of the first second-dimension column, the overall ²D selectivity was tuned. The practical advantages of this tunable system were evaluated by measuring the ²D relative retention of 60 target analytes and by focusing on critical pairs of compounds in a commercial roasted coffee as a real-world sample. The analysis posed some challenges related to ²D chromatographic resolution of critical pairs, for example: 2-methyl-3-hydroxy-4-pyrone (maltol)/ 1-methylpyrrole-2-carboxaldehyde and 4-hydroxy-2,5-dimethyl-3(2H)-furanone /2-acetylpyrrole.

Another example of how hyphenation can improve the informative potential of GC×GC was presented by Tranchida et al. [46], who combined high-performance liquid chromatography (LC) and GC×GC with fast quadrupole MS (qMS) in order to characterize cold-pressed sweet orange oil and bergamot essential oils. Pre-separation was performed by means of a LC×GC system with a 100 mm × 3 mm ID × 5 μm *d_p* silica column operated under gradient elution with hexane/MTBE as the mobile phase at 0.35 mL/min constant flow. Fractions were collected on the basis of their polarity: hydrocarbons were collected from 1.5 to 3 min (525 μL); sweet orange oil oxygenated compounds were collected from 7.3 to 14 min (2345 μL); and bergamot oil oxygenated compounds were collected from 7.5 to 13 min (1925 μL). Prior to GC×GC-MS analysis, fractions were reduced to 100 μL under a gentle stream of nitrogen.

Experimental results for the sweet orange oil were straightforward, as the authors stated in their concluding remarks: a total number of 219 analytes were identified, compared to 50 solutes assigned by using 1D-GC-MS as reference method. Of the identified analytes, 169 reported a spectrum similarity match probability greater than 90% and a difference in Linear Retention Index ($\pm I^T$) of 5 units or less. In addition, 38 analytes had not been reported previously. A total of 195 analytes were identified in bergamot oil, compared with 64 assigned by 1D-GC-MS. Of the identified analytes, 171 reported a spectrum similarity match probability greater than 90% and a $\pm I^T$ of 5 or less. Twenty new compounds were tentatively identified and were shown to be present in bergamot oil for the first time.

322

323 **Hyphenation with sample preparation**

324 In reviewing *-omics* investigations of food sensory quality, sample preparation deserves a dedicated
325 section, as this is one of the bottlenecks of the entire analytical process. In order to deliver a consistent and
326 meaningful picture of volatiles and semi-volatiles, including sensory-active analytes, a sample preparation
327 technique must provide: (a) *ad-hoc* tuning of the extraction selectivity by modifying physico-chemical
328 characteristics of the extractants and sampling conditions (time, temperature, and volume/mass of the
329 extraction phase); (b) flexibility in terms of extraction efficiency/capability, because the absolute amount
330 extracted directly affects method performance in terms of Limit of Determination (LOD) and Limit of
331 Quantitation (LOQ); (c) extraction methods based on mild interactions to limit artifact formation, thus

sorption (i.e., partition) should be preferred vs. adsorption as the extraction mechanism; and (d) the possibility of full integration and automation of the extraction process, thus including sample preparation as an additional dimension in the analytical platform [47-49].

In this context, well established extraction procedures, such as Solvent Assisted Flavor Evaporation (SAFE), Simultaneous Distillation Extraction (SDE), Cold Finger distillation (CF), hydro distillation, Solid Phase Extraction (SPE), and Supercritical fluids Extraction (SFE), which have been used for many years by flavor chemists, have been replaced, whenever possible, by automated approaches, because these techniques have limited possibilities for hyphenation with the analytical platform.

Above all, headspace extraction approaches have regained strong interest because of demonstrated capabilities on a wider range of applications in the food field. These techniques, also classified as High Concentration Capacity Headspace Techniques (HCC-HS) [50], offer an elective route for satisfactory throughput headspace sampling. They are based on either a static or dynamic accumulation of volatiles on polymers operating in sorption and/or adsorption. Selectivity and extraction capability can be tuned *ad hoc* to meet the requirements for a given application, by selecting appropriate polymers, their physical state, and volume. In particular, HS-SPME and Headspace Sorptive Extraction (HSSE) are the most widely-used static HCC-HS approaches, are easy to standardize, and can be integrated in the separation system. Dynamic headspace sampling (D-HS) can be considered as a valid alternative, being able to increase sensitivity and achieve higher concentration factors [47], although a careful tuning of sampling parameters is necessary to avoid breakthrough and to obtain a representative picture of volatiles without discriminations [51-56].

Headspace Solid Phase Microextraction is undoubtedly the most popular of the HCC-HS techniques and its hyphenation with GC×GC platforms is well documented in a number of applications, some listed in **Table 1**. Rochat et al. [57] investigated sulfur containing odorants of beef by extracting volatiles directly from the oven headspace while a piece of meat was roasted. The application required the sensitivity of GC×GC-TOFMS coupled to an enrichment technique in the extraction step, because sulfur compounds are potent odorants that often occur at trace levels. Volatiles from vapors were extracted by inserting an SPME silicone fiber (PDMS - 100 µm) for 10 min inside a glass condenser installed on the down-stream of an *ad-hoc* designed tubular ventilated oven. An additional extraction, aimed at enriching trace and sub-trace analytes, was conducted with an organomercurial derivative of the N-hydroxysuccinimide activated agarose gel for affinity chromatography (Affi-gel 501, BioRad Reinach, Switzerland). The stationary phase enabled the selective isolation of mercaptans (SH) that were successively eluted in different fractions and also assayed by panelists. Fractions exhibiting the most intense odor were mixed and submitted to HS-SPME sampling (PDMS - 100µm) before GC×GC-TOFMS analysis. This approach enabled identification of seven impact odorants from among 69 sulfur derivatives (23 thiophenes, 19 thiazoles, and 27 mercaptans, sulfide, and isothiocyanate derivatives), of which six exhibited the highest impact in the roast-beef top note: 2-

367 methyl-3-mercapto-1-propanol was characterized by *beef broth, meaty, onion juice* notes; 3-
368 (methylthio)thiophene by *alliaceous, sulfurous, rubbery, gassy, coffee*; (±)-2-Methyl-3-[(2-
369 methylbutyl)thio]furan by *meaty, green, weak, sulfurol, meaty*; 2-phenylthiophene by *vague, rubbery,*
370 *weak*; 3-phenylthiophene by *meaty, rubbery*; 4-isopropylbenzenethiol by *mushroom, alliaceous, cardboard*;
371 and 4-(methylthio)benzenethiol by a *rubbery, weak* note. With the exception of 2-methyl-3-mercapto-1-
372 propanol, which also was reported to occur in wine, the other compounds were identified for the first time
373 in beef and none of them had been previously mentioned in survey listing food aroma compounds from
374 TNO[58].

375 A study by Chin et al. [43] discusses the advantages, in terms of detection limit for GC-O screening,
376 obtained by using cumulative HS-SPME as sample preparation for wine aroma assessment. Such an
377 experimental design presents challenging aspects: the difficulty of automation and, from the GC-O
378 perspective, of performing replicate assays or dilution experiments. The proposed method included twelve
379 contemporary samplings with two different fiber coatings, followed by successive GC injections delayed
380 over time.

381 In a study focused on hazelnut aroma characterization, Nicolotti et al. [59] moved a step forward and
382 proposed a quantitative method based on Multiple Headspace Extraction (MHE) with SPME. The approach,
383 whose advantages will be discussed in more detail in the application section, provided information not only
384 on the concentrations of analytes, but showed interesting fingerprinting potential because only minimal
385 differences were detectable in the chemical pattern when the headspace linearity condition is matched
386 [59]. Thanks to the high sensitivity of the GC×GC-MS, the number of matched peaks within 2D
387 chromatograms only decreased from 100% with the 1.500 g sample to 73% with the 0.100 g sample. More
388 precisely, 73 unknown and 17 known analytes were lost by sampling 0.100 g and only a few odor-active
389 compounds and one key aroma compound (i.e., 2-acetyl-1-pyrroline) fell below the method LOD. MHE-
390 SPME-GC×GC-MS applied to food-end products prepared with hazelnut paste (*Gianduja* paste: sugar,
391 vegetable oil, hazelnuts, cocoa, nonfat milk, vanilla flavorings) provided also a measure of the actual
392 release of some key-odorants (2,3-pentanedione, 5-methyl-(*E*)-2-hepten-4-one, (*E*)-2-octenal, 2,5-dimethyl-
393 3-ethyl-pyrazine, 2,6-dimethyl-3-ethyl-pyrazine, phenylacetaldehyde, (*E*)-2-decenal, 3-methylbutanoic acid,
394 2-phenylethanol, and acetylpyrrole).

395 Gogus et al. [60] investigated the effect of roasting time on the volatiles of *Pistacia terebinthus* L. fruit,
396 growing wild in Turkey. Whole fruits were pan roasted and successively submitted to direct thermal
397 desorption (DTD) followed by GC×GC-TOFMS analysis. DTD, although of interest because of the ease of use
398 and the possibility of automation, in this specific application presents some drawbacks related to the
399 thermal exposure of the matrix during desorption time. Non-volatile constituents undergo thermal
400 degradation producing a pattern of volatile derivatives that interfere with the univocal identification of
401 those formed exclusively during the pan roasting.

402 Villire et al. [61] investigated the potential of SPME (applied as HS or as in-solution sampling), Dynamic
403 Headspace extraction (D-HS) with polar adsorbents (i.e., Tenax), and Purge and Trap (P&T) to provide
404 representative extracts of French cider for GC-O screenings. The HS-SPME fiber coating polymers in
405 particular were investigated. Car/PDMS was found to be the most suitable coating to obtain representative
406 HS profiles of cider odor. Experimental designs for fiber selection and extraction conditions (time and
407 temperature) were oriented by the similarity score and representativeness of the chromatographic profile
408 combined with a sensory assay conducted by 12 panelists, who were asked to evaluate the gaseous phase
409 trapped on a glass syringe. Aromagrams obtained by GC-O revealed 36 and 24 odorant zones for the two
410 cider samples, which were subsequently investigated by GC×GC-TOFMS.

411 Van der Wat et al. [62] adopted PDMS traps, i.e., Multi Channel Silicone Rubber Traps (MCT), to
412 characterize the volatile fraction of rosemary (*Rosmarinus officinalis* L.) from two different geographical
413 origins, Tunisia and South Africa.

414 A study by Cordero et al. [63] on a volatile fraction isolated from dried milk reported a systematic
415 investigation on the effectiveness of different and complementary hyphenated and automated sampling
416 techniques, based on either sorption and adsorption, or a combination of them, with the aim of
417 qualitatively and quantitatively screening volatiles and semi-volatiles of dry milk powders, especially
418 focusing on sensory-active analytical targets (key-aroma compounds and off-odorants). Approaches
419 investigated, most of them carried out automatically, were SPME, Stir Bar Sorptive Extraction (SBSE) and
420 Head Space Sorptive Extraction (HSSE) with silicone and dual phase extraction media, and Dynamic
421 Headspace sampling (D-HS) with silicone sorbents or polar adsorbents like Tenax TA™. The information for
422 analytes extracted by headspace and in-solution sampling were compared to evaluate whether a given
423 orthogonal approach was advantageous to describe the sensory properties of the investigated samples. The
424 sample matrix investigated, i.e., dry milk powders (whole and nonfat milk), posed some challenges because
425 of the wide range of volatility (Vapor pressures - Vp), polarity (LogP values from 0.3 to 8), water solubility,
426 and concentration of the most significant analytes, which required that both powders and reconstituted
427 liquids would be analyzed for a reliable characterization of the final aroma profile. **Figure 3** reports the 2D
428 patterns of a whole dried milk sample and its linear saturated aldehydes from C-6 to C-18 obtained by
429 Dynamic Headspace sampling with PDMS traps, i.e. D-HS-PDMS, as well as the 2D pattern of lactones
430 resulting from a Headspace Sorptive Extraction sampling with PDMS stir-bar, i.e. HSSE-PDMS. 2D plots are
431 obtained by selecting diagnostic *m/z* fragments (i.e., 57,82,95 *m/z* for aldehydes and 55,71,99 *m/z* for
432 lactones) from the Total Ion Current (TIC) (top of image) by scripting with CLIC™ Expression (GC-Image, LLC
433 Lincoln NE, USA) on the software platform [63].

434 Among the investigated techniques, HSSE and SBSE were highly effective for sensomics because of their
435 high concentration factors, allowing them to provide highly descriptive profiles as well as analyte amounts
436 suitable for GC-O screenings, even with high odor threshold (OT) markers or potent odorants in sub-trace

amounts. Therefore, the approach represents a possible bridge between classic extraction procedures (LLE, SDE and SAFE) and more popular approaches such as SPME.

It should, however, be stressed that for an exhaustive and truly comprehensive characterization of key-aroma compounds, classical procedures of isolation of the volatile fraction performed on suitably high amounts of sample matrix may be required. Kiefl et al. [64] introduced a useful parameter to evaluate the performance of an analytical method to measure concentrations at the odor threshold level by considering the LOQ. The parameter, defined as Limit of Odor Activity Value (LOAV), corresponds to the ratio between analyte odor threshold (OT) and method LOQ. By definition, a LOAV greater than 1 indicates a sensitive method that gives an effective and quantitative odorant assessment above the odor threshold, whereas a LOAV less than 1 indicates the concentration limit under which an odorant can be identified but not accurately quantified.

Bi-dimensional data elaboration challenges

Comprehensive two-dimensional chromatography offers unequaled information on compositional characteristics of complex samples, but the data size and complexity make data analysis to extract information a challenging problem. Cross-sample analysis in this specific field of application aims, for example: (a) to classify samples on the basis of their sensory profile; (b) to obtain chemical fingerprints to correlate sample characteristics with those of reference samples; (c) to monitor progressive or cyclical changes as a function of a specific technological/enzymatic treatment; (d) to cluster similar samples; and/or (e) to discover informative markers of botanical/geographical origin.

The most relevant *features* (i.e., analytical entities characterized by detector or mass spectral intensities) for a particular cross-sample analysis sometimes are related to trace analytes and/or unidentified compounds. Thus, a productive investigation strategy should start with a non-targeted approach to extract and analyze all information that may be relevant. However, non-targeted analysis requires dedicated software(s) and skillful analysts to perform chemometrics procedures to reduce and rationalize data processing outputs. On the other hand, an extended untargeted processing would be unnecessary for those applications where, for instance, a bio-guided assay (e.g. GC-O) preliminarily targets/tags specific retention regions as meaningful to describe the sensory properties of a sample.

Most of the studies reviewed here have adopted targeted approaches, by first identifying analytes on the basis of their EI-MS fragmentation pattern and relative retention (by linear retention indices - I^T s) and successively comparing relative distributions across samples. Multivariate Analysis (MVA) is frequently adopted in post-processing, with both unsupervised and supervised approaches to select those variables within a set that better “describe” the problem under investigation.

Vaz-Freire et al. [65] investigated the effects of two extraction methods used in the production of Extra Virgin Olive (EVO) oils (i.e., metal hammer-decanter vs. traditional metal hammer-press line) on the aroma

472 compounds from Portuguese varieties *Galega*, *Carrasquenha*, and *Cobrançosa*. Bi-dimensional patterns
473 obtained by HS-SPME sampling and GC×GC-TOFMS from freshly extracted oils were processed by a *region*
474 *feature* approach performed with open-source software (ImageJ™, Wayne Rasband, National Institute of
475 Health, USA). *Region features* consist of datapoint clusters in the chromatographic plane (e.g., summing the
476 intensities at all datapoints in each region) that characterize meaningful chromatographic structures. In this
477 application, the authors covered the entire chromatographic space with rectangles of equal size (1000 s in
478 the ¹D and 2 s in the ²D) in which analytes are present. The response from each rectangle was collected and
479 used for cross-sample analysis. ANOVA after Tukey validation confirmed the consistency of the region
480 feature results, in terms of cumulative response, when compared to 2D peak distributions. Principal
481 Component Analysis (PCA) was able to cluster samples according to their botanical origin and to locate the
482 most informative regions where discriminating analytes eluted.

483 Schmarr and Bernhardt [66] analyzed volatile patterns, including some aroma-active compounds, from
484 apple, pear, and quince fruits and adopted an advanced profiling analysis approach for cross-samples
485 comparison. Volatiles, sampled by HS-SPME, were successively analyzed with GC×GC-qMS to generate a
486 unique informative data matrix for each single analysis. Data was converted to a jpeg image by open-source
487 software (ImageJ™, Wayne Rasband, National Institute of Health, USA) and processed with a *peak-region*
488 *feature* approach commonly adopted for 2D gel electrophoresis. This approach consisted of a sequence of
489 pre-processing operations (images were aligned and summed) that produced a single chromatogram
490 representative of all of the constituents in all samples. **Figure 4** summarizes the workflow of the proposed
491 method. The boundaries that delineated each peak were recorded as a region in a template. The template
492 was geometrically mapped back to each chromatogram and detector responses (intensities) extracted and
493 compared across the samples set. Feature matching was performed by retention-times mapping; MS data
494 were not included as a matching restriction. Post-processing and data interpretation was by Hierarchical
495 Cluster Analysis (HCA) and PCA on the *peak-region features*. The different fruits formed clear clusters, and
496 sub-clusters were formed by pear and some apple varieties.

497 The same approach also was adopted to differentiate microoxygenation (MOX) treatments and varietal and
498 technological effects on Pinot noir, Cabernet Sauvignon, and Dornfelder wines of the 2007 vintage [67]. The
499 authors identified peak-regions that could be used to discriminate between the different MOX treatments
500 and the loadings of individual aroma compounds suggested a set of markers for the MOX-induced
501 modifications of volatiles.

502 Smart Templates™ with peak-region features were developed by Reichenbach and co-workers [68] and
503 used to characterize the volatile fraction of coffee and *Juniper* samples. After preprocessing, including peak
504 detection, peaks that could be matched reliably across all chromatograms were identified. These reliable
505 peaks, with mass spectral matching rules, were used to build a registration template, which was then used
506 to determine the geometric transforms to align each of the chromatograms. After alignment, the

507 chromatograms were summed to create a composite chromatogram. In three chromatograms of coffee
508 samples, about 1700 peaks were detected, about half of which were reliable. They manually drew a mesh
509 of about 1100 regions which were combined with the registration peaks to create a feature template that
510 could be matched to individual chromatograms. The geometry of the reliable peak matching was used to
511 transform the regions in order to maintain their positions relative to the reliable peaks. The features were
512 sifted by intensity, standard deviation, and relative standard deviation to select relevant features, but MVA
513 was not used because of the small number of samples. Many of the indicated compounds were known
514 botanical, technological, and/or aromatic markers for coffee. For the analysis of the five chromatograms of
515 juniper samples, there were about 100 reliable peaks and 727 peak-regions were drawn.

516 Bordiga et al. [69] developed a pixel-based approach on 2D raw data from HS-SPME-GC×GC-TOFMS analysis
517 of volatiles from different Muscat wines from Piedmont stored at different temperatures for six months (5,
518 15 and 25°C). The method, classified as a *pointwise* approach, enabled point-by-point (or in imaging terms
519 pixel-by-pixel) chromatographic comparisons; each datapoint was used as a feature and the datapoint
520 features at the same retention times were implicitly matched.

521 Cordero et al. [70] investigated the volatile fraction of roasted hazelnuts from different botanical and
522 geographical origins with HS-SPME-GC×GC-qMS and non-targeted cross-comparisons based on *peak*
523 *features*, with Comprehensive Template Matching (CTM) fingerprinting. Templates for peak matching were
524 obtained with two different approaches. In the first approach, they aligned and summed the
525 chromatograms and then created a feature template with the 411 peaks detected in the cumulative
526 chromatogram. This template was matched to each individual chromatogram, with peak-matching rates
527 ranging from 68% to 79%. In the second approach, they performed a sequential template matching that
528 used both retention-time patterns and mass spectral matching criteria. At each matching step, unmatched
529 peaks were added to the comprehensive template. At the end of the sequence, the comprehensive
530 template was matched to each chromatogram and all peaks matching at least two chromatograms were
531 retained in a *consensus template*. The consensus template contained 422 peaks and the matching rates
532 ranged from 52% to 78%, with 196 peaks matching for all nine chromatograms. For both peak matching
533 methods, the feature fingerprints of samples from nine geographic regions were sifted for the largest
534 normalized intensities and many of the indicated compounds were known to have a role in determining
535 sensory properties.

536 In a successive study, Kiefl et al. [71] validated the CTM fingerprinting approach on a series of hazelnut
537 samples from different origins and technological treatments and concluded that an appropriate setting of
538 data elaboration parameters (peak detection thresholds based on SNR, retention-time search windows, MS
539 match factor thresholds, and template thresholds) would limit false positive/negative matching and
540 improve the reliability of non-targeted cross-comparison of samples. The validated method successfully
541 elucidated the generation of volatile compounds during roasting in a set of 23 hazelnut samples, in which

11 roasting markers were identified. The results showed that the release of key aroma compounds produced specific profiles as a function of variety/origin of hazelnut samples. Purcaro et al. [72] adopted CTM fingerprinting followed by supervised MVA to identify the blueprint of regulated defects of Extra Virgin Olive (EVO) oils. Nineteen olive oil samples, including 5 reference standards obtained from the International Olive Oil Council and 14 commercial samples were submitted for sensory evaluation by a panel, prior to an analysis in two laboratories using different instrumentation, column sets, and software elaboration packages in view of a cross-validation of the entire methodology. A first classification of samples, based on untargeted *peak features*, was obtained on raw data from two different column combinations (apolar × polar and polar × apolar) by PCA. However, to improve effectiveness and specificity of the classification, peak features were reliably identified (261 compounds) on the basis of the MS spectrum and linear retention index matching, then subjected to successive pair-wise comparisons based on 2D patterns, which revealed peculiar distributions of chemicals correlated with samples' sensory classification. The most informative compounds were identified and collected in a *blueprint* of specific defects (or combination of defects) successively adopted to discriminate Extra Virgin from defected oils (i.e. *lampante* oil) with the aid of a supervised approach, i.e., Partial Least Squares-Discriminant Analysis (PLS-DA). In the last step, the principle of Sensomics, assigning higher information potential to analytes with lower odor threshold, proved to be successful, and a much more powerful discrimination of samples was obtained in view of a sensory quality assessment.

560

561 **Applications**

Up to now, the characterization of key odorants, such as off-odorants and character impact compounds, has been one of the most important applications of gas chromatography in food aroma analysis. A few hundred aroma compounds have been identified by using GC-O in more than 100 different foods [12]. The challenge of detecting trace amounts of highly active odorants in complex food matrices has been a motivating force for the development of more sensitive methods with higher peak capacity and with increased linear detector response. Therefore, with the first commercially available GC×GC instruments, researchers transferred and developed new analytical methods to characterize key odorants [19,30,63,64,73,74]. Accordingly, this review considers first the characterization of key odorants using GC×GC (e.g., detailed profiling of key-odorants); second, the differentiation of aromas by correlation of key odorant fingerprints with sensory data; and finally, the identification of marker compounds to predict aroma profiles. Investigations on essential oils of interest in food applications are not discussed here.

573

574 *Key-odorants profiling*

The two potent aroma compounds methional and sotolon were identified and quantified at 35 µg/kg and 85 µg/kg in sour cream and dairy spread extracts obtained by SAFE and CF distillation using GC×GC-TOFMS

577 in two-dimensional and one-dimensional mode [31]. The comparison of one and two-dimensional
578 separations showed that coelution of these components could be minimized and sensitivity improved and,
579 moreover, the elution order of homologous series of aroma compounds was a valuable tool for the
580 identification of unknowns.

581 Wine was one of the first food matrices to be investigated for the odorant 3-isobutyl-2-methoxypyrazine
582 (IBMP) with GC×GC. Ryan et al. [73] quantified IBMP in Sauvignon blanc wine by HS-SPME-GC×GC-TOFMS
583 and using the ²H₃-isotopomer as internal standard. A limit of detection of 1.95 ng/L, similar to already
584 existing methods using 1D-GC, was reported; however, comparably less time for sample treatment was
585 needed. By investigating the same analyte in Sauvignon blanc wine, Schmarr et al. [74] concluded that
586 GC×GC separation alone might not be enough for proper chromatographic resolution and that more
587 powerful mass spectrometers such as (high-resolution)-TOFMS compared to quadrupole MS are needed for
588 additional mass spectrometric resolution.

589 Fresh lemon juice and thermally stressed lemon-flavored beverages were analyzed with GC×GC-FID to
590 identify the citral degradation products *p*-cymen-8-ol and *p*-methylacetophenone, which play a significant
591 role in off-flavor development [75]. Identification of just 24 volatile compounds could be achieved by co-
592 chromatography revealing that the lack of further structural information, such as mass spectral data, makes
593 identification tedious and time consuming [75].

594 Poisson et al. [76] quantified 3-methyl-2-butene-1-thiol in coffee brew, an odorant which may play a key
595 role for the overall aroma of freshly ground coffee, by using HS-SPME-GC×GC-TOFMS. On average, 0.12
596 µg/L were determined, while saving time compared to SPE combined with H/C MDGC-MS analysis, thanks
597 to the enhanced peak capacity of GC×GC, which minimized co-elution and increased SNR.

598 Forty-seven odorants with an OAV > 1 were identified in Chardonnay wine by external calibration in model
599 wine using HS-SPME-GC×GC-TOFMS [77]. Compared to other studies analyzing Chardonnay wine odorants
600 with GC-MS, a higher number of compounds with OAV > 1 were found, thereby demonstrating the
601 capability of GC×GC-TOFMS to profile odorants more effectively.

602

603 *Correlation of odorant fingerprints with sensory data*

604 The identification of food odorants by additionally employing GC-Olfactometry to characterize smell often
605 has been used in combination with GC×GC-MS analysis [57,61,78-80]. For example, Rochat et al. [57]
606 identified more than 25 odor-active sulfur compounds with a SNIF value > 50 % in roasted beef and Villire
607 et al. [61] more than 20 odorants in French ciders by combining GC-O with GC×GC-TOFMS. Breme et al. [81]
608 located 44 odorants in an extract of Indian cress using GC-O and the vocabulary-intensity-duration of
609 elementary odors (VIDEO) sniff technique to identify 22 of them with HS-SPME-GC×GC-TOFMS, including
610 (E)-hex-2-enal (*fruity*) and diethyl trisulfide (*alliaceous, sulfury, cabbage*), which were found to have the
611 highest odor impact. Thirty odor-active compounds have been identified in cereal coffee brew via GC-O and

612 Aroma Extract Dilution Analysis (AEDA) by Majcher et al. [79], 17 of them with an OAV > 1 after
613 quantitation with Stable Isotope Dilutions Assays (SIDA) and Standard Addition (SA) with GC×GC-TOFMS.
614 GC-O, AEDA, and SIDA/SA combined with GC×GC-TOFMS have also been used to decode the aroma of
615 fermented and fried soy tempeh formed by the specific ratios of 2-acetyl-1-pyrroline, 2,6-dimethyl-3-ethyl-
616 pyrazine, dimethyltrisulfide, methional, 2-methylpropanal, and (E,E)-2,4-decadienal [80]. Finally,
617 Maikhunthod and Marriott [44] identified limonene, 1,8-cineole, terpinen-4-ol, estragole, and *trans*-
618 anethole as main aroma compounds in dried fennel seeds by GC-O/NIF and HS-SPME-GC×GC-TOFMS.

619 Although hundreds of volatiles might be detected with GC×GC-MS, GC-O guides the attention to a
620 few key odorants and, hence, makes olfactometry still an essential tool [78]. The combination of GC-
621 O/GC×GC-MS, however, requires the correlation of retention times, respective to retention indices, to
622 define a small retention-time window on the 2D chromatographic plane where the potential aroma
623 compound elutes. Even if the same columns and chromatographic parameters are used, such retention-
624 time windows could be large enough to present too many peaks for an unambiguous mass spectral
625 identification. Rochat et al. [78], for example, correlated I_s^T of GC-MS-O, MDGC-O, and GC×GC-MS for the
626 identification of 23 shrimp aroma compounds with a NIF value above 50 %.

627 Chin et al. [41] detected more than 200 volatile compounds with SPE extraction and GC×GC-TOFMS in
628 brewed coffee, Merlot, and a white wine blend, and located 19 odor-active chromatographic zones with a
629 SNIF value > 50 % for brewed coffee and 14 for Merlot by using GC-O×GC-FID. The odor-active compound
630 for each chromatographic zone, however, could not always be determined, because the odor descriptors of
631 the analytes obtained from literature and online databases eluting in these zones did not match with the
632 descriptors of the GC-O×GC-FID experiment. Eyres et al. [82] located with GC-O Charm analysis between 38
633 and 71 odor-active zones in the spicy fraction of four different hops and could identify the corresponding
634 aroma compound in just 13 of 25 zones investigated, leaving 12 zones as unknown.

635 These studies show that correlating the odor perceived at the GC sniffing port from one-
636 dimensional separation with the mass spectrometric data from two-dimensional separation is challenging,
637 because more than one of the peaks spread along the second dimension may fit the recorded odor quality
638 or a peak with recorded odor quality does not match any known compound. For this reason, Eyres et al.
639 [82] employed MDGC-O to locate odorants by sniffing in the second dimension as well. Especially if
640 unknown odorants have to be identified, sniffing in the second dimension is mandatory. D'Acampora
641 Zellner et al. [83] hyphenated the second dimension with a sniffing port (GC×GC-O), however, peaks elute
642 within the millisecond range and, although the modulation period might be lengthened, a high breathing
643 rate would be needed to sufficiently resolve the peaks for detection by the human nose. For this reason,
644 this technique is not yet established. GC(O)×GC-MS is used to bypass the correlation of retention times in
645 the first dimension [39,61]. Thus, the number of peaks can be constrained by setting smaller retention-time

646 windows for identification while inter-instrumental variations are excluded, but this still might not provide
647 sufficient confirmatory evidence to assign peak identity in the second dimension.

648 Although GC×GC-MS is considered to be a complementary tool for the characterization of key odorants by
649 providing enhanced peak capacity and sensitivity to facilitate the identification of trace amounts of
650 odorants coeluting with highly abundant odor-inactive compounds [39,78,82], H/C MDGC is still the
651 method of choice for the unambiguous identification of unknown odorants in the second dimension. As
652 discussed previously, H/C MDGC can be combined with GC×GC to cut a modulation sequence rather than a
653 conventional retention-time window, thus giving the possibility to identify the unknown odorant within the
654 same run on a second dimension.

655 Beyond characterizing a few key odorants, an increasing number of publications aim at profiling the entire
656 set of volatiles of a food with GC×GC and correlate the data with sensory analysis (**Figure 5**). For this
657 reason, the sample preparation and GC×GC analysis could be optimized to assess quantitatively (or at least
658 relatively) the concentrations or area ratios and to correlate the Chemical Odor Code to sensory data such
659 as Quantitative Descriptive Analysis (QDA) [84] or projective mapping [85]. Multivariate statistical methods
660 such as PCA, Discriminant Analysis (DA), Artificial Neural Network (ANN), and Multi Dimensional Scaling
661 (MDS) or simple calculative operations can be used to establish correlations and to develop models to
662 predict the results of sensory tests with a glimpse of the Chemical Odor Code. Validation of such models
663 (**Figure 5**), e.g., by mixing model solutions and performing sensory evaluation, is essential because there is
664 no other test system than our olfactory sense which can closely mimic human odor perception [86].

665 The sensory characteristics of Cabernet Sauvignon wines from different locations were studied by
666 correlating the distribution of over 350 volatiles with QDA data from 16 aroma attributes [84]. A trained
667 panel with 18 assessors provided the sensory data and the volatiles were analyzed by HS-SPME-GC×GC-
668 TOFMS. Wines characterized as *fruity* and *vegetal herbaceous* could be well differentiated by correlating
669 the *fruity* note (among others) to δ -octalactone, vitispirane, γ -decalactone, and γ -octalactone, and the
670 *vegetal/herbaceous* note to IBMP, which smells like *bell pepper*. The model suggested eucalyptol and
671 hydroxy-citronellol as being important for the *eucalypt* and *mint* aroma attributes; furan and benzene
672 derivatives were positively correlated with the aroma perception of *oak*; and the *floral* characteristic was
673 connected with dihydro- α -ionone and sesquiterpenes like α - and β -calacorene. Although these correlations
674 sound reasonable, no confirmatory evidence, e.g., by spiking experiments, was given to show cause-effect
675 relations between sensory attributes and proposed compounds [84]. The same research group used the
676 identical approach to investigate the role of yeast, canopy, and site on the composition and sensory
677 characteristics of Western Australian Cabernet Sauvignon wines [87].

678 Inui et al. [86] brewed beer with five different aroma hops and studied the aroma compounds by
679 correlating QDA data of the 6 attributes, *floral*, *herbal*, *citric*, *spicy*, *ester*, and *sylvan*, from 5 trained
680 panelists, with GC×GC-TOFMS analytical data. For example, the hop “Tradition” showed a high score in the

681 ester character, “Perle” was high in *sylvan* character, and “Cascade” beer showed the highest intensity,
682 especially in *floral*, *citric*, and *spicy* notes. Multivariate analysis with PCA indicated the correlation of 67
683 compounds from 297 volatiles detected with the above six sensory descriptors. However, the authors
684 suggested running further experiments to prove the results by sensory experiments [86].

685 Purcaro et al. [72] correlated the peak fingerprint of 19 different olive oil samples with the sensory
686 properties classified in *musty*, *vinegary*, *fusty*, *mold*, *rancid*, and *fruity*. The number of more than 400
687 volatiles was reduced and the normalized 2D peak volume of statistically significant peaks were submitted
688 first to PCA then to PLS-DA also using the ratio of normalized 2D peak volume and odor threshold. It was
689 shown that better classification and hence correlation with the sensory attributes was obtained when the
690 odor-activity value was considered by including odor-thresholds.

691 Fifteen different samples from three hazelnut cultivars of different geographic origin and roasting
692 degree were analyzed by SAFE-GC×GC-TOFMS and Stable Isotope Dilution Analysis (SIDA) in order to profile
693 over 20 odorants quantitatively [85]. These analytical data were correlated with sensory data from a
694 projective mapping experiment with 20 panelists visualizing aroma differences and similarities on a 2D
695 plane. The resulting aroma map was matched with the Chemical Odor Code by simple calculative
696 operations: odorants exceeding the threshold concentration were first selected by calculating OAVs ($OAV \geq$
697 1), then these odorants were grouped according to their aroma attributes assuming synergistic effects, and
698 finally concentrations were iteratively drawn on a x/y coordinate system to find the pattern with the
699 highest aroma map similarity. The model suggested that the *roasty*, *nutty* aroma of optimally roasted
700 hazelnuts was developed if both 5-methyl-(E)-2-hepten-4-one and 3-methyl-4-heptanone exceeded 450
701 µg/kg, whereas the sum of 2-acetyl-/2-propionyl-1-pyrroline, 2,5(6)-dimethyl-3-ethylpyrazine, and 2,3-
702 diethyl-5-methylpyrazine should not to exceeded 400 µg/kg to avoid an *over-roasted* smell. The hypothesis
703 was successfully tested by mixing the proposed odorants in deodorized sunflower oil and submitting these
704 model mixtures again to projective mapping (**Figure 6**). In **Figure 6A**, the results of sensory analysis are
705 shown and compared to the sensory evaluation of the model mixtures obtained by correlation (**Figure 6B**).
706 Three main clusters could be defined: raw hazelnuts on the left; a group of optimally roasted samples with
707 a nutty, roasty smell in the upper part; and the over-roasted samples on the bottom right corner. Further
708 sensory experiments to substantiate the model by studying odorant interactions on the basis of odorants’
709 natural concentrations were conducted and provided deep insights into the mechanisms of the aroma
710 development in hazelnuts [86].

711 These examples show that correlating analytical data with sensory analysis is challenging because
712 the mechanisms of odor perception driven by interactions of key odorants is more complex than single
713 statistical methods can delineate. Obviously, no standard approach is available for this purpose because the
714 statistical methods and experimental designs used in these studies are unique for each subject. Hence,
715 GC×GC can provide a more detailed picture of the Chemical Odor Code compared to 1D-GC, but to

understand how this Code is translated into an aroma profile, the developed models should be validated by studying sensory effects in model solutions. In this view, GC×GC based methods (a) should be optimized to characterize the whole set of key odorants as quickly as possible so the sample preparation or aroma extraction should be soft but representative, and (b) should provide the best separation for selected volatiles and, if necessary, detect traces of them within a wide linear range. Furthermore, quantitative profiling would facilitate the validation of developed models by re-engineering respective aromas. Nicolotti et al. [59] recently proposed a quantitation method based on Multiple Headspace Extraction (MHE) with SPME for the quantitative profiling of 19 relevant odorants and technological markers of the roasting process of hazelnuts. This quantitation method should be useful and provide additional insights into the release of volatiles from the food matrix [59].

726

727 *Aroma profile prediction*

Comprehensive two-dimensional gas chromatography also can be successfully used to profile the volatile fractions of a food to identify marker compounds that may be related to the sensory properties by statistical means. Commonly, no cause-effect relations between marker compounds and the aroma are established and, hence, such compounds are used as indicators for quality assessment.

The effect of microoxygenation on the volatile fraction of three different red wines was studied by HS-SPME-GC×GC-qMS, and PCA statistical analysis [67]. Six alcohols, 2 aldehydes, 1 ketone, 8 esters, β-damascenone, and 1,1,6-trimethyl-1,2-dihydronaphthalene were proposed as markers of microoxygenation, which is supposed to improve the flavor of red wine.

The volatile pattern of Albion and Juliette strawberries were compared by HS-SPME-GC×GC-TOFMS analysis [88]. In Albion, γ-decalactone, methyl butanoate, methyl hexanoate, (E)-hex-2-enal, and (E)-nerolidol were the most abundant volatile constituents; in Juliette, the most abundant volatile constituents were (E)-hex-2-enal, (E)-nerolidol (14.6%), mesifuran, (E)-hex-2-enyl acetate, and linalool. The detection of 2,5-dimethyl-4-hydroxy-(2H)-furan-3-one only in Juliette and the higher area percentage of 2,5-dimethyl-4-methoxy-(2H)-furan-3-one were considered to be correlated to the enhanced sweetness of Juliette strawberries.

The concept of profiling the volatile fraction and identifying aroma marker compounds also was applied to fresh picked and stored strawberries [89], Chinese liquor Maotai [90], different basil (*Ocimum basilicum* L.) cultivars grown under conventional and organic conditions [91], different honeys [92,93], malolactically fermented Trincadeira wine [94], roasted hazelnuts [70], roasted *Pistacia terebinthus* L. fruit [60], different South African Pinotage wines [95], Muscat wines [69], and green, oolong and black teas [96].

747

748 **Conclusions and future perspectives**

For many analyses, GC×GC is the technique currently offering the highest peak capacity, but its potential in many fields has not yet been fully exploited. In the flavor field, this technique has been shown to be

751 valuable for highly detailed characterization of food volatiles; to study the composition of complex volatile
752 fractions, very often consisting of hundreds of components (e.g., coffee, tea or cocoa); to detect key-
753 odorants and explain their formation from precursors; and to understand the interaction/relationship with
754 flavor perception, personal behavior, and health. Although not always indispensable, because the relevance
755 of information cannot be reduced to the number of components that can be separated, GC×GC offers
756 higher separation power and sensitivity that can be fundamental for (a) accurate aroma fingerprints of
757 complex samples (e.g., processed food) to be correlated with sensory perceptions and, as a consequence,
758 sensory qualities, and (b) a better aroma blueprints of food, i.e., the distribution of key-aroma compounds,
759 in particular when present in trace amounts. GC×GC is especially promising for flavor research on particular
760 problems not solvable with conventional techniques, but it still requires a degree of sophistication that is
761 rather high, which adversely affects its routine use.

762 Recent instrumental advances have integrated olfactometry (GC(O)×GC-MS) and H/C MDGC in
763 conventional GC×GC-MS platforms. Although Olfactometry does not increase dramatically the complexity
764 of the system, simultaneous GC×GC/HCMDGC may require a level of sophistication such that is adoptable
765 only in highly specialized laboratories. Flow modulation is a relatively simple and cheap technology and,
766 once some of its technological problems are overcome, e.g., with the introduction of reverse flow
767 modulation systems or flexible loop-capillary [24,36,38], can contribute to adoption of GC×GC for routine
768 analysis.

769 Crucial to GC×GC advancement is its relationship with other fundamental steps, i.e., sample preparation,
770 analyte isolation, and data elaboration. The present trend, in general, and in particular with GC×GC, is to
771 achieve, where possible, a full integration with sample preparation in order to include it as a further
772 dimension of a fully automatic separation platform. A strong effort is, therefore, underway to combine
773 sample preparation and GC×GC on-line. In the field of food volatile fraction and aroma characterization,
774 this trend has contributed greatly to a renewed interest in headspace sampling, in all its modes (static (S-
775 HS), dynamic (D-HS) and high concentration capacity techniques (HCC-HS)), because it can easily be
776 integrated on-line with the analytical instrumentation. In addition, the possibility to adopt concentration
777 materials operating on different principles (sorption and adsorption) and of different chemical natures with
778 D-HS and HCC-HS can be very useful as a preliminary selective step when specific classes of compounds
779 have to be analyzed. The HS sampling success is confirmed by the fact that about 80% of the articles quoted
780 in this review adopt HS-SPME, not only because of its undeniable effectiveness, but also because of its ease
781 of integration in a total analysis system. This trend also has resulted in the development of fully automatic
782 and versatile purge and trap and D-HS sampling systems operating in series with different accumulation
783 phases, thus extending the applicability of GC×GC and GC(O)×GC to samples where high concentration
784 factors are required.

785 Data elaboration is the step of the analytical process that is expected to be the object of the most radical
786 evolution in the next few years. Two important main trends in this respect are:
787 1) Improved data pre-processing will reduce ambiguities in 2D peak detection and peak area/volume
788 determination with better standardized and more widely accepted algorithms (as is the case for 1D GC).
789 2) Data analysis will extract more information. The latter trend merits a more detailed comment with
790 respect to applications in flavor research. At present, conventional data elaboration is mainly based on
791 targeted profiling, which is limiting because it excludes all other data on known and unknown components
792 deriving from the comprehensive separation. This approach can be satisfactory for those applications
793 where retention regions of targets that are significantly representative of the sensory properties of a
794 sample are known, e.g., by GC-O. On the other hand, a truly effective elaboration strategy implies the
795 adoption of a non-targeted approach (advanced fingerprinting), in which information useful to characterize
796 the investigated aroma, and as a consequence food, can be extracted from all data made available from the
797 chromatographic separation. However, at present, advanced fingerprinting requires dedicated software(s)
798 and external chemometrics procedures to reduce and rationalize data processing outputs. As is
799 concurrently developing in metabolomics, further processing tools in this direction are expected to become
800 more effective and mainstream.

801

802 References

- 803 1. Oliver SG, Winson MK, Kell DB, Baganz F (1998) Systematic functional analysis of the yeast
804 genome. Trends in Biotechnology 16: 373-378
- 805 2. de Vos CHR, Tikunov Y, Bovy A G, Hall RD (2008) Flavour metabolomics: Holistic versus targeted
806 approaches in Flavour research. In Expression of Multidisciplinary Flavour Science. Blank I, Wust M,
807 Yeretzian C Proceedings of the 12th Weurman Symposium, Zürcher Hochschule für Angewandte and
808 Institut Für Chemie und Biologische Chemie, Interlaken, Switzerland, 573-580
- 809 3. Herrero M, Simó C, García-Cañas V, Ibáñez E, Cifuentes A (2012) Foodomics: MS-based strategies in
810 modern food science and nutrition. Mass Spectrometry Reviews 31: 49-69
- 811 4. Schieberle P, Hofmann T (2011). Mapping the Combinatorial Code of Food Flavors by Means of
812 Molecular Sensory Science Approach. Food Flavors, CRC Press: 413-438
- 813 5. Charve J, Chen C, Hegeman AD, Reineccius GA, (2011) Evaluation of instrumental methods for the
814 untargeted analysis of chemical stimuli of orange juice flavor. Flavour and Fragrance Journal 26 (6)
815 429-440
- 816 6. Harrigan G, Goodacre R (2003) Metabolic profiling: its role in biomarker discovery and gene
817 function analysis. Kluwer Academic Publishers, Boston, 2003
- 818 7. Grosch W (2001) Evaluation of the Key Odorants of Foods by Dilution Experiments, Aroma Models
819 and Omission. Chem Senses 26: 533-545
- 820 8. Cordero C, Rubiolo P, Bicchi C, (2008) Group-type and fingerprint analysis of roasted food matrices
821 (coffee and hazelnut samples) by comprehensive two-dimensional gas chromatography. Journal of
822 Agricultural and Food Chemistry. 56 (17): 7655-7666
- 823 9. Firestein S.(2001) How the olfactory system makes sense of scents. Nature, 413: 211-218

10. Fleischer J, Breer H, Strotmann J (2009) Mammalian olfactory receptors. *Frontiers in Cellular Neuroscience* 3: 1-10
11. Breer H, Fleischer J, Strotmann J (2006) The sense of smell: multiple olfactory subsystems. *Cellular and Molecular Life Sciences* 63: 1465–1475
12. Dunkel A, Steinhaus M, Kotthoff M, Nowak B, Krautwurst D, Schieberle P, Hofmann T (2014) Nature's Chemical Signatures in Human Olfaction: A Foodborne Perspective for Future Biotechnology. *Angewandte Chemie International Edition*. doi: 10.1002/anie.201309508
13. Teranishi R, Wick EL, Hornstein IEds. (1999) *Flavor Chemistry- Thirty Years of Progress*, Kluwer Academic Plenum Publishers, New York
14. Cortes H, Ed. (1990) *Multidimensional Chromatography: Techniques and Applications*. Dekker, New York
15. Mondello L, Lewis A C, Bartle KD Eds. (2001) *Multidimensional Chromatography*. Wiley, Chichester, UK
16. Marriott P J, Kinghorn RM (2001) Multidimensional and comprehensive multidimensional gas chromatography. In *Gas Chromatographic Techniques and Applications*. Handley A J, Adlard E R, Eds. Sheffield Academic Press, Sheffield, UK 260-297.
17. Delahunty C M, Eyres G, Dufour JP (2006) Gas chromatography olfactometry. *Journal of Separation Science* 29:2107-2125
18. Marriott PJ, Eyres GT, Dufour JP (2009) Emerging opportunities for flavor analysis through hyphenated gas chromatography. *Journal of Agricultural and Food Chemistry*. 57(21):9962-9971
19. Marriott PJ, Chin ST, Maikhunthod B, Schmarr HG, Bieri S (2012) Multidimensional gas chromatography. *Trends in Analytical Chemistry* 34:1-21
20. Chambers E, Koppel K (2013) Associations of Volatile Compounds with Sensory Aroma and Flavor: The Complex Nature of Flavor. *Molecules* 18 (5):4887-4905
21. Harynuk T, Górecki T, de Zeeuw J (2005) Overloading of the second-dimension column in comprehensive two-dimensional gas chromatography. *Journal of Chromatography A*. 1071: 21-27
22. Giddings J C (1995) Sample dimensionality: a predictor of order-disorder in component peak distribution in multidimensional separation. *Journal of Chromatography A*. 703: 3-15.
23. Seeley JV, Seeley SK, (2013) Multidimensional gas chromatography: Fundamental advances and new applications. *Analytical Chemistry* 85(2): 557-578
24. Tranchida PQ, Sciarrone D, Dugo P, Mondello L (2012) Heart-cutting multidimensional gas chromatography: A review of recent evolution, applications, and future prospects. *Analytica Chimica Acta*. 716: 66-75
25. Bertsch W. (1999) Two-dimensional gas chromatography. Concepts, instrumentation, and applications - Part 1: fundamentals, conventional two-dimensional gas chromatography, selected applications. *Journal of High Resolution Chromatography*. 22: 647-665.
26. Liu Z, Phillips JB (1991) Comprehensive two-dimensional gas chromatography using an on-column thermal modulator interface. *Journal of Chromatographic Science*. 29(6): 227-231
27. Phillips JB, Xu J (1995) Comprehensive multi-dimensional gas chromatography. *Journal of Chromatography A*. 703: 327-334
28. Kidwell DA, Riggs LA (2004) Comparing two analytical methods: minimal standards in forensic toxicology derived from information theory. *Forensic science international*. 145: 85-96.
29. Davis JM, Giddings JC (1983) Statistical-theory of component overlap in multicomponent chromatograms. *Analytical Chemistry*. 55: 418–424
30. Davis JM, Samuel C (2000) The need for two-dimensional gas chromatography: extent of overlap in one-dimensional gas chromatograms. *Journal of High Resolution Chromatography*. 23: 235-244.

31. Adahchour M, Beens J, Vreuls RJJ, Batenburg AM, Rosing EAE, Brinkman UATH (2002) Application of solid-phase micro-extraction and comprehensive two-dimensional gas chromatography (GC x GC) for flavour analysis. *Chromatographia*. 55: 361-367
32. Adahchour M, Van Stee LLP, Beens J, Vreuls RJJ, Batenburg MA, Brinkman UATH (2003) Comprehensive two-dimensional gas chromatography with time-of-flight mass spectrometric detection for the trace analysis of flavour compounds in food. *Journal of Chromatography A*. 1019: 157-172
33. Engel W, Bahr W, Schieberle P (1999) Solvent Assisted Flavour Evaporation - a New Versatile Technique for the Careful and Direct Isolation of Aroma Compounds from Complex Food Matrices. *European Food Research and Technology*. 209: 237-241
34. Widder S, Sen A, Grosch W, (1991) Changes in the flavour of butter oil during storage. *Zeitschrift für Lebensmittel-Untersuchung und Forschung*. 193(1): 32-35
35. Tranchida PQ, Purcaro G, Dugo P, Mondello L (2011) Modulators for comprehensive two-dimensional gas chromatography. *Trends in Analytical Chemistry*. 30(9): 1437-1461
36. Manzano P, Diego JC, Bernal JL, Nozal MJ, Bernal J (2014) Comprehensive two-dimensional gas chromatography coupled with static headspace sampling to analyze volatile compounds: Application to almonds. *Journal of Separation Science*. 37(6): 675-683
37. Tranchida PQ, Purcaro G, Visco A, Conte L, Dugo P, Dawes P, Mondello L (2011) A flexible loop-type flow modulator for comprehensive two-dimensional gas chromatography. *Journal of Chromatography A*. 1218: 3140-3145
38. Tranchida PQ, Franchina FA, Zoccali M, Bonaccorsi I, Cacciola F, Mondello L (2013) A direct sensitivity comparison between flow-modulated comprehensive 2D and 1D GC in untargeted and targeted MS-based experiments. *Journal of Separation Science*. 36(17): 2746-2752
39. Chin ST, Eyres GT, Marriott PJ (2012) System design for integrated comprehensive and multidimensional gas chromatography with mass spectrometry and olfactometry. *Analytical Chemistry*. 84(21): 9154-9162
40. Maikhunthod B, Morrison PD, Small DM, Marriott PJ (2010) Development of a switchable multidimensional/comprehensive two-dimensional gas chromatographic analytical system. *Journal of Chromatography A*. 1217: 1522-1529
41. Chin ST, Eyres GT, Marriott PJ (2011) Identification of potent odourants in wine and brewed coffee using gas chromatography-olfactometry and comprehensive two-dimensional gas chromatography. *Journal of Chromatography A*. 1218: 7487-7498
42. Pollien P, Ott A, Montigon F, Baumgartner M, Munoz-Box R, Chaintreau A (1997) Hyphenated Headspace-Gas Chromatography-Sniffing Technique: Screening of Impact Odorants and Quantitative Aromagram Comparisons. *Journal of Agricultural and Food Chemistry*. 45: 2630-2637
43. Chin ST, Eyres GT, Marriott PJ (2012) Cumulative solid phase microextraction sampling for gas chromatography-olfactometry of Shiraz wine. *Journal of Chromatography A*. 1255: 221-227
44. Maikhunthod B, Marriott PJ (2013) Aroma-impact compounds in dried spice as a quality index using solid phase microextraction with olfactometry and comprehensive two-dimensional gas chromatography. *Food Chemistry*. 141(4): 4324-4332
45. Mommers J, Pluimakers G, Knooren J, Dutriez T, van der Wal S(2013) Tunable secondary dimension selectivity in comprehensive two-dimensional gas chromatography. *Journal of Chromatography A*. 1297: 179-185
46. Tranchida PQ, Zoccali M, Bonaccorsi I, Dugo P, Mondello L, Dugo G(2013) The off-line combination of high performance liquid chromatography and comprehensive two-dimensional gas

- chromatography-mass spectrometry: A powerful approach for highly detailed essential oil analysis. *Journal of Chromatography A*. 1305: 276-284
47. Bicchi C, Cordero C, Liberto E, Sgorbini B, Rubiolo P (2012) Headspace sampling in flavor and fragrance field. *Comprehensive Sampling and Sample Preparation*. 4: 1-25
48. Lord HL, Pfannkoch EA (2012) Sample preparation automation for GC injection. *Comprehensive Sampling and Sample Preparation*. 2: 597-612
49. Risticevic S., Vuckovic D., Lord H.L., Pawliszyn J., Solid-phase microextraction, 2012 ,*Comprehensive Sampling and Sample Preparation*, 2, 419-460
50. Bicchi C, Cordero C, Liberto E, Sgorbini B, Rubiolo P (2008) Headspace sampling of the volatile fraction of vegetable matrices. *Journal of Chromatography A*. 1184: 220-233
51. Sgorbini B, Cagliero C, Cordero C, Liberto E, Rubiolo P, Ruosi MR, Bicchi C (2012) New medium-to-high polarity twister coatings for liquid and vapour phase sorptive extraction of matrices of vegetable origin. *Journal of Chromatography A*. 1265: 39-45
52. Snow NH, Slack GC (2002) Head-space analysis in modern gas chromatography. *Trends in Analytical Chemistry*. 21:608-617.
53. Baltussen E, Cramers CA ,Sandra P (2002) Sorptive sample preparation - a review. *Analytical and Bioanalytical Chemistry*. 373: 3-22.
54. Chary NS, Fernandez-Alba AR (2012) Determination of volatile organic compounds in drinking and environmental waters. *Trends in Analytical Chemistry*. 32: 60-75.
55. Tholl D, Boland W, Hansel A, Loreto, Röse USR, Schnitzler (2006) Practical approaches to plant volatile analysis. *The Plant Journal*. 45: 54-560
56. Augusto F, Leite E, Lopes A, Zini CA (2003) Sampling and sample preparation for analysis of aromas and fragrances. *Trends in Analytical Chemistry*. 22: 160-169
57. Rochat S, de Saint Laumer JY, Chaintreau A (2007) Analysis of sulfur compounds from the in-oven roast beef aroma by comprehensive two-dimensional gas chromatography. *Journal of Chromatography A*. 1147: 85-94
58. VCF00: database of volatile compounds in food. Boelens aroma chemical information service, BACIS, Zeist, 1999.
59. Nicolotti L, Cordero C, Cagliero C, Liberto E, Sgorbini B, Rubiolo P, Bicchi C (2013) Quantitative fingerprinting by headspace-Two-dimensional comprehensive gas chromatography-mass spectrometry of solid matrices: Some challenging aspects of the exhaustive assessment of food volatiles. *Analytica Chimica Acta*. 798: 115-125
60. Gogus F, Ozel MZ, Kocak D, Hamilton JF, Lewis AC (2011) Analysis of roasted and unroasted *Pistacia terebinthus* volatiles using direct thermal desorption-GC×GC-TOF/MS. *Food Chemistry*. 129(3): 1258-1264
61. Villire A, Arvisenet G, Lethuaut L, Prost C, Serot T (2012) Selection of a representative extraction method for the analysis of odourant volatile composition of French cider by GC-MS-O and GC × GC-TOF-MS. *Food Chemistry*. 131(4): 1561-1568
62. Van Der Wat L, Dovey M, Naude Y, Forbes PBC (2013) Investigation into the aroma of rosemary using multi-channel silicone rubber traps, off-line olfactometry and comprehensive two-dimensional gas chromatography-mass spectrometry. *South African Journal of Chemistry*. 66: 21-26
63. Cordero C, Cagliero C, Liberto E, Nicolotti L, Rubiolo P, Sgorbini B, Bicchi C (2013) High concentration capacity sample preparation techniques to improve the informative potential of two-dimensional comprehensive gas chromatography-mass spectrometry: Application to sensomics. *Journal of Chromatography A*. 1318,: 1-11

64. Kiefl J, Pollner G, Schieberle P (2013) Sensomics analysis of key hazelnut odorants (*Corylus avellana* L. 'Tonda Gentile') using comprehensive two-dimensional gas chromatography in combination with time-of-flight mass spectrometry (GC×GC-TOF-MS). *Journal of Agricultural and Food Chemistry*. 61(22): 5226-5235
65. Vaz-Freire LT, da Silva MDRG, Freitas AMC (2009) Comprehensive two-dimensional gas chromatography for fingerprint pattern recognition in olive oils produced by two different techniques in Portuguese olive varieties Galega Vulgar, Cobrançosa e Carrasquenha. *Analytica Chimica Acta*. 633(2): 263-270
66. Schmarr HG, Bernhardt J (2010) Profiling analysis of volatile compounds from fruits using comprehensive two-dimensional gas chromatography and image processing techniques. *Journal of Chromatography A*. 1217: 565-574
67. Schmarr HG, Bernhardt J, Fischer U, Stephan A, Muller P, Durner D (2010) Two-dimensional gas chromatographic profiling as a tool for a rapid screening of the changes in volatile composition occurring due to microoxygenation of red wines. *Analytica Chimica Acta*. 672:114-123
68. Cordero C, Liberto E, Bicchi C, Rubiolo P, Reichenbach SE, Tian X, Tao Q (2010) Targeted and non-targeted approaches for complex natural sample profiling by GC×GC-qMS. *Journal of Chromatographic Science*. 48(4): 251-261
69. Bordiga M, Rinaldi M, Locatelli M, Piana G, Travaglia F, Coisson JD, Arlorio M (2013) Characterization of Muscat wines aroma evolution using comprehensive gas chromatography followed by a post-analytic approach to 2D contour plots comparison. *Food Chemistry*. 140: 57-67
70. Cordero C, Liberto E, Bicchi C, Rubiolo P, Schieberle P, Reichenbach SE, Tao Q (2010) Profiling food volatiles by comprehensive two-dimensional gas chromatography coupled with mass spectrometry: Advanced fingerprinting approaches for comparative analysis of the volatile fraction of roasted hazelnuts (*Corylus avellana* L.) from different origins. *Journal of Chromatography A*. 1217: 5848-5858
71. Kiefl J, Cordero C, Nicolotti L, Schieberle P, Reichenbach SE, Bicchi C (2012) Performance evaluation of non-targeted peak-based cross-sample analysis for comprehensive two-dimensional gas chromatography-mass spectrometry data and application to processed hazelnut profiling. *Journal of Chromatography A*. 1243: 81-90
72. Purcaro G, Cordero C, Liberto E, Bicchi C, Conte LS (2014) Toward a definition of blueprint of virgin olive oil by comprehensive two-dimensional gas chromatography. *Journal of Chromatography A*. 1334: 101-111
73. Ryan D, Watkins P, Smith J, Allen M, Marriott P (2005) Analysis of methoxypyrazines in wine using headspace solid phase microextraction with isotope dilution and comprehensive two-dimensional gas chromatography. *Journal of Separation Science*. 28 (9-10): 1075-1082
74. Schmarr H-G, Ganß S, Koschinski S, Fischer U, Riehle C, Kinnart J, Potouridis T, Kutyrev M (2010) Pitfalls encountered during quantitative determination of 3-alkyl-2-methoxypyrazines in grape must and wine using gas chromatography - mass spectrometry with stable isotope dilution analysis. Comprehensive two-dimensional gas chromatography - mass spectrometry and on-line liquid chromatography-multidimensional gas chromatography - mass spectrometry as potential loopholes. *Journal of Chromatography A*. 1217: 6769-6777
75. Komura H (2006) Comprehensive two-dimensional gas chromatographic analysis of commercial lemon-flavored beverages. *Journal of Separation Science*. 29: 2350-2356
76. Poisson L, Hug C, Baggenstoss J, Blank I, Kerler (2011) Emerging analytical techniques for the assessment of aroma relevant sulfur compounds in coffee. *ACS Symposium Series 1068*: 77-92

77. Welke JE, Zanusi M, Lazzarotto M, Alcaraz Zini C (2014) Quantitative analysis of headspace volatile compounds using comprehensive two-dimensional gas chromatography and their contribution to the aroma of Chardonnay wine. *Food Research International* 59: 85-99
78. Rochat S, Egger J, Chaintreau A (2009) Strategy for the identification of key odorants: Application to shrimp aroma. *Journal of Chromatography A*. 1216: 6424-6432
79. Majcher MA, Klensporf-Pawlik D, Dziadas M, Jelen HH (2013) Identification of Aroma Active Compounds of Cereal Coffee Brew and Its Roasted Ingredients. *Journal of Agricultural and Food Chemistry* .61:2648-2654.
80. Jelen H, Majcher M, Ginja A, Kuligowski M (2013) Determination of compounds responsible for tempeh aroma. *Food Chemistry*.141(1): 459-465
81. Breme K, Tournayre P, Fernandez X, Meierhenrich UJ, Brevard H, Joulain D, Berdague JL (2009) Identification of Odor Impact Compounds of *Tagetes minuta* L. Essential Oil: Comparison of Two GC-Olfactometry Methods. *Journal of Agricultural and Food Chemistry*. 57: 8572-8580
82. Eyres GT, Marriott PJ, Dufour J-P (2007) Comparison of Odor-Active Compounds in the Spicy Fraction of Hop (*Humulus lupulus* L.) Essential Oil from Four Different Varieties. *Journal of Agricultural and Food Chemistry*. 55: 6252-6261
83. d'Acampora Zellner B, Casilli A, Dugo P, Dugo G, Mondello L (2007) Odor fingerprint acquisition by means of comprehensive two-dimensional gas chromatography-olfactometry and comprehensive two-dimensional gas chromatography/mass spectrometry. *Journal of Chromatography A*. 1141: 279-286
84. Robinson AL, Boss PK, Heymann H, Solomon PS, Trengove RD (2011) Development of a sensitive non-targeted method for characterizing the wine volatile profile using headspace solid-phase microextraction comprehensive two-dimensional gas chromatography time-of-flight mass spectrometry. *Journal of Chromatography A*. 1218: 504-517
85. Kiefl J, Schieberle P (2013) Evaluation of process parameters governing the aroma generation in three hazelnut cultivars (*Corylus avellana* L.) by correlating quantitative key odorant profiling with sensory evaluation. *Journal of Agricultural and Food Chemistry*. 61: 5236-5244
86. Inui T, Tsuchiya F, Ishimaru M, Oka K, Komura H (2013) Different beers with different hops. Relevant compounds for their aroma characteristics. *Journal of Agricultural and Food Chemistry*. 61: 4758-4764
87. Robinson AL, Boss PK, Heymann H, Solomon PS, Trengove RD (2011) Influence of Yeast Strain, Canopy Management, and Site on the Volatile Composition and Sensory Attributes of Cabernet Sauvignon Wines from Western Australia. *Journal of Agricultural and Food Chemistry*. 59: 3273-3284
88. Samykanno K, Pang E, Marriott PJ (2013) Chemical characterisation of two Australian-grown strawberry varieties by using comprehensive two-dimensional gas chromatography-mass spectrometry. *Food Chemistry*. 141(3): 1997-2005
89. Williams A, Ryan D, Olarte Guasca A, Marriott PJ, Pang E (2005) Analysis of strawberry volatiles using comprehensive two-dimensional gas chromatography with headspace solid-phase microextraction. *Journal of Chromatography B: Analytical Technologies in the Biomedical and Life Sciences*. 817: 97-107
90. Zhu S, Lu X, Ji K, Guo K, Li Y, Wu C, Xu G (2007) Characterization of flavor compounds in Chinese liquor Moutai by comprehensive two-dimensional gas chromatography/time-of-flight mass spectrometry. *Analytica Chimica Acta*. 597 (2):340-348.

91. Klimankova E., Holadova K., Hajslova J., Cajka T., Poustka J., Koudela M (2008) Aroma profiles of five basil (*Ocimum basilicum* L.) cultivars grown under conventional and organic conditions. Food Chemistry. 107: 464-472
92. Cajka T, Hajslova J, Pudil F, Ridelova K (2009) Traceability of honey origin based on volatiles pattern processing by artificial neural networks., Journal of Chromatography A. 1216: 1458-1462
93. Stanimirova I, Ustun B, Cajka T, Ridelova K, Hajslova J, Buydens LMC, Walczak B (2010) Tracing the geographical origin of honeys based on volatile compounds profiles assessment using pattern recognition techniques. Food Chemistry. 118: 171-176
94. Pietra Torres M, Cabrita MJ, Gomes Da Silva MDR, Palma V, Costa Freitas AM (2011) The impact of malolactic fermentation on the volatile composition of the trincadeira wine variety. Journal of Food Biochemistry. 35(3): 898-913
95. Weldegergis BT, Villiers AD, McNeish C, Seethapathy S, Mostafa A, Gorecki T, Crouch AM (2011) Characterisation of volatile components of Pinotage wines using comprehensive two-dimensional gas chromatography coupled to time-of-flight mass spectrometry (GC x GC-TOFMS). Food Chemistry. 129(1): 188-199
96. Zhang L, Zeng Z, Zhao C, Kong H, Lu X, Xu G (2013) A comparative study of volatile components in green, oolong and black teas by using comprehensive two-dimensional gas chromatography-time-of-flight mass spectrometry and multivariate data analysis. Journal of Chromatography A. 1313: 245-252
97. Mondello L, Casilli A, Tranchida PQ, Dugo P, Costa R, Festa S, Dugo G (2004) Comprehensive multidimensional GC for the characterization of roasted coffee beans. Journal of Separation Science. 27: 442-450
98. Ryan D, Shellie R, Tranchida P, Casilli A, Mondello L, Marriott P (2004) Analysis of roasted coffee bean volatiles by using comprehensive two-dimensional gas chromatography-time-of-flight mass spectrometry. Journal of Chromatography A. 1054: 57 -65
99. Adahchour M, Brandt M, Baier H-U, Vreuls RJJ, Batenburg AM, Brinkman UATH (2005) Comprehensive two-dimensional gas chromatography coupled to a rapid-scanning quadrupole mass spectrometer: Principles and applications. Journal of Chromatography A. 1067: 245-254
100. Eyres G, Dufour JP, Hallifax G, Sotheeswaran S, Marriott PJ (2005) Identification of character-impact odorants in coriander and wild coriander leaves using gas chromatography-olfactometry (GCO) and comprehensive two-dimensional gas chromatography-time-of-flight mass spectrometry (GC x GC-TOFMS). Journal of Separation Science. 28: 1061-1074
101. Mondello L, Costa R, Tranchida PQ, Dugo P, Presti ML, Festa S, Fazio A, Dugo G (2005) Reliable characterization of coffee bean aroma profiles by automated headspace solid phase microextraction-gas chromatography-mass spectrometry with the support of a dual-filter mass spectra library. Journal of Separation Science. 28: 1101-1109
102. Williams A, Ryan D, Olarte Guasca A, Marriott PJ, Pang E (2005) Analysis of strawberry volatiles using comprehensive two-dimensional gas chromatography with headspace solid-phase microextraction. Journal of Chromatography B: Analytical Technologies in the Biomedical and Life Sciences. 817: 97-107
103. Cardeal ZL, Gomes Da Silva MDR, Marriott PJ (2006) Comprehensive two-dimensional gas chromatography/mass spectrometric analysis of pepper volatiles. Rapid Communications in Mass Spectrometry. 20(19): 2823-2836
104. Chaintreau A, Rochat S, de Saint Laumer J-Y (2006) Re-investigation of sulfur impact odorants in roast beef using comprehensive two-dimensional GC-TOF-MS and the GC-SNIF technique. Developments in Food Science. 43: 601-604

105. de Saint Laumer J-Y, Chaintreau A (2006) Prediction of gas-chromatographic retention indices as a tool for identification of sulfur odorants. *Developments in Food Science*. 43: 597-600
106. Bianchi F, Careri M, Conti C, Musci M, Vreuls R (2007) Comparison of comprehensive two-dimensional gas chromatography-time-of-flight mass spectrometry and gas chromatography-mass spectrometry for the qualitative characterisation of roasted barley by solid-phase microextraction. *Journal of Separation Science*. 30(4): 527-533
107. Cajka T, Hajslova J, Cochran J, Holadova K, Klimankova E (2007) Solid phase microextraction-comprehensive two-dimensional gas chromatography-time-of-flight mass spectrometry for the analysis of honey volatiles. *Journal of Separation Science*. 30(4): 534-546
108. Eyres G, Marriott PJ, Dufour J-P (2007) The combination of gas chromatography-olfactometry and multidimensional gas chromatography for the characterisation of essential oils. *Journal of Chromatography A*. 1150: 70-77
109. Rocha SM, Coelho E, Zrostlikova J, Delgadillo I, Coimbra MA (2007) Comprehensive two-dimensional gas chromatography with time-of-flight mass spectrometry of monoterpenoids as a powerful tool for grape origin traceability. *Journal of Chromatography A*. 1161: 292-299
110. Cardeal ZL, de Souza PP, Silva MDRGd, Marriott PJ (2008) Comprehensive two-dimensional gas chromatography for fingerprint pattern recognition in cachaça production . *Talanta*. 74(4): 793-799
111. Klimankova E, Holadova K, Hajslova J, Cajka T, Poustka J, Koudela M (2008) Aroma profiles of five basil (*Ocimum basilicum* L.) cultivars grown under conventional and organic conditions. *Food Chemistry*. 107(1): 464-472
112. Cardeal ZL, Marriott PJ (2009) Comprehensive two-dimensional gas chromatography-mass spectrometry analysis and comparison of volatile organic compounds in Brazilian cachaça and selected spirits. *Food Chemistry*. 112(3): 747-755
113. de Souza PP, Cardeal ZL, Augusti R, Morrison P, Marriott PJ (2009) Determination of volatile compounds in Brazilian distilled cachaça by using comprehensive two-dimensional gas chromatography and effects of production pathways. *Journal of Chromatography A*. 1216: 2881-2890
114. Humston EM, Zhang Y, Brabeck GF, McShea A, Synovec RE (2009) Development of a GC×GC-TOFMS method using SPME to determine volatile compounds in cacao beans. *Journal of Separation Science*. 32(13): 2289-2295
115. Lojzova L, Ridelova K, Hajslova J, Zrostlikova J, Schurek J, Cajka T (2009) Alternative GC-MS approaches in the analysis of substituted pyrazines and other volatile aromatic compounds formed during Maillard reaction in potato chips. *Analytica Chimica Acta*. 641: 101-109
116. Humston EM, Knowles JD, McShea A, Synovec RE (2010) Quantitative assessment of moisture damage for cacao bean quality using two-dimensional gas chromatography combined with time-of-flight mass spectrometry and chemometrics. *Journal of Chromatography A*. 1217: 1963-1970
117. Silva I, Rocha SM, Coimbra MA, Marriott PJ (2010) Headspace solid-phase microextraction combined with comprehensive two-dimensional gas chromatography time-of-flight mass spectrometry for the determination of volatile compounds from marine salt. *Journal of Chromatography A*. 1217: 5511-5521
118. Stanimirova I, Ustun B, Cajka T, Ridelova K, Hajslova J, Buydens LMC, Walczak B (2010) Tracing the geographical origin of honeys based on volatile compounds profiles assessment using pattern recognition techniques. *Food Chemistry*. 118(1): 171-176

119. Chin S-T, Eyres GT, Marriott PJ (2011) Identification of potent odourants in wine and brewed coffee using gas chromatography-olfactometry and comprehensive two-dimensional gas chromatography. *Journal of Chromatography A*. 1218: 7487-7498
120. Weldegergis BT, Villiers AD, McNeish C, Seethapathy S, Mostafa A, Gorecki T, Crouch AM (2011) Characterisation of volatile components of Pinotage wines using comprehensive two-dimensional gas chromatography coupled to time-of-flight mass spectrometry (GC×GC-TOFMS). *Food Chemistry*. 129(1): 188-199
121. Omar J, Alonso I, Olivares M, Vallejo A, Etxebarria N (2012) Optimization of comprehensive two-dimensional gas-chromatography (GC×GC) mass spectrometry for the determination of essential oils. *Talanta*. 88: 145-151
122. Langos D, Granvogl M, Schieberle P (2013) Characterization of the key aroma compounds in two Bavarian wheat beers by means of the sensomics approach. *Journal of Agricultural and Food Chemistry*. 61(47): 11303-11311
123. Rivellino SR, Hantao LW, Risticvic S, Carasek E, Pawliszyn J, Augusto F (2013) Detection of extraction artifacts in the analysis of honey volatiles using comprehensive two-dimensional gas chromatography. *Food Chemistry*. 141(3): 1828-1833
124. Willner B, Granvogl M, Schieberle P (2013) Characterization of the key aroma compounds in bartlett pear brandies by means of the sensomics concept. *Journal of Agricultural and Food Chemistry*. 61(40): 9583-9593
125. Bordiga M, Piana G, Coisson JD, Travaglia F, Arlorio M (2014) Headspace solid-phase micro extraction coupled to comprehensive two-dimensional with time-of-flight mass spectrometry applied to the evaluation of Nebbiolo-based wine volatile aroma during ageing. *International Journal of Food Science and Technology*. 49(3): 787-796
126. Dugo G, Franchina FA, Scandinaro MR, Bonaccorsi I, Cicero N, Tranchida PQ, Mondello L (2014) Elucidation of the volatile composition of marsala wines by using comprehensive two-dimensional gas chromatography. *Food Chemistry*. 142: 262-268

1169 **Caption to Tables**

1170 **Table 1:** Overview of GC×GC applications in food sensory quality characterization

1171

1172 **Caption to Figures**

1173 **Figure 1:** Schematic diagram of the switchable targeted MDGC/GC GC×GC system. DS: Deans switch; CT:
1174 cryotrap; 1D: first dimension column; 2DS: short second dimension column (for GC×GC mode) terminated
1175 at Flame Ionization detector FID 1; 2DL: long second dimension column (for targeted MDGC mode)
1176 terminated at Flame Ionization detector FID 2 [40].

1177

1178 **Figure 2:** Instrument schematic of the integrated GC × GC/MDGC system with olfactory and mass spectral
1179 detection. DS: Deans switch; ES: effluent splitter; AUX: auxiliary pressure port; 1D: first dimension column;
1180 2DL: second dimension long column; 2DS: second dimension short column; CT: cryotrap; SSI: split/splitless
1181 injector; FID: flame ionization detector [39].

1182

1183 **Figure 3:** 2D pattern of a whole dried milk sample submitted to D-HS with PDMS packing and HSSE-PDMS
1184 sampling. (a) The TIC trace of the sample headspace (D-HS-PDMS), (b) the SIM trace of linear saturated
1185 aldehydes (57,82,95 m/z) and (c) SIM trace for lactones (55,71,99 m/z) recovered by HSSE-PDMS. SIM
1186 images were obtained by scripting with CLICTM Expression (GC-Image, LLC Lincoln NE, USA) [63].

1187

1188 **Figure 4:** (1) Samples have been prepared and analyzed by HS-SPME-GC×GC-qMS; (2) 2D GC
1189 chromatograms have been transformed into 32-bit images; (3) 2D GC images were stored in Delta2DTM
1190 software; (4) Positional correction (warp vectors) resulted in image congruency (dual channel overlay color
1191 code: blue = image1, orange = image2 and black = overlap); (5) Volatiles map as a result of project-wide 2D
1192 GC image fusion; (6) Detected spot consensus; (7) Spot consensus boundaries were applied to all 2D GC
1193 images for gray level integration; (8) Gray level integration results in quantitative data which can be
1194 summarized in volatile profiles (blue – low amount, black – average amount, orange – large amount of
1195 volatile) [66].

1196

1197 **Figure 5:** schematic procedure of the correlation of volatile compound fingerprints with sensory data
1198 obtained from the same samples using “omics” techniques [86].

1199


1200 **Figure 6:** (A): Consensus perceptual map of the raw (0) as well as 12, 23 and 30 min roasted and industrially
1201 roasted (no number) hazelnut cultivars Akçakoca (A), ‘Gentile’ (G), and ‘Romana’ (R). Duplicate samples are
1202 presented as controls (c). (B): Consensus perceptual map of respective odorant model mixtures suggested
1203 by correlation based on sunflower oil containing between four to eight odorants [85].

Table 1

Year	Authors	Instrumental platform(s)	Sample preparation	Data elaboration	Food matrix	Ref
2002	Adahchour et al.	GC×GC-FID	HS-SPME	Targeted profiling	flavour analysis	31
2003	Adahchour et al.	GC×GC-TOFMS; GC-TOFMS	SAFE and Cold Finger (CF) distillation	Targeted profiling	dairy products	32
2004	Mondello et al.	GC×GC-qMS	HS-SPME	Targeted profiling	roasted coffee	97
2004	Ryan et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	roasted coffee	98
2005	Adahchour et al.	GC×GC-qMS	Direct injection	Targeted profiling	flavour analysis	99
2005	Eyres et al.	GC-O; GC×GC-TOFMS	EO hydrodistillation	Targeted profiling	<i>Coriandrum sativum</i> ; <i>Eryngium foetidum</i>	100
2005	Mondello et al.	GC×GC-qMS	HS-SPME	Targeted profiling	roasted coffee	101
2005	Ryan et al.	GC×GC-NPD; GC×GC-TOFMS	HS-SPME	Targeted profiling; Quantitative fingerprinting	wine	73
2005	Williams et al.	ES-GC×GC-FID	HS-SPME	Targeted profiling	strawberry	102
2006	Cardeal et al.	GC×GC-TOFMS; GC×GC-qMS	HS-SPME	Targeted profiling	pepper and peppercorn	103
2006	Chaintreau et al.	GC-O; GC×GC-TOFMS	HS-SPME, affinity chromatography	Targeted profiling	roast-beef	104
2006	de Saint Laumer & Chaintreau	GC×GC-TOFMS	HS-SPME, affinity chromatography	Targeted profiling	sulfur odorants	105
2006	Komura	GC×GC-FID	LLE	Targeted profiling	lemon flavored beverages	75
2007	Bianchi et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	roasted barely	106
2007	Cajka et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	honey	107
2007	d'Acampora Zellner et al.	GC×GC-O/qMS	Direct injection	Targeted profiling	commercial perfumes	83
2007	Eyres et al.	GC-O; GC×GC-TOFMS	EO hydrodistillation	Targeted profiling	<i>Coriandrum sativum</i> and <i>Humulus lupulus</i> essential oils	108
2007	Eyres et al.	GC-O; GC×GC-TOFMS	EO hydrodistillation	Targeted profiling	hop essential oils	82
2007	Rocha et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	grape	109
2007	Rochat et al.	GC-O; GC×GC-TOFMS	HS-SPME, affinity chromatography	Targeted profiling	roast-beef	57
2008	Cardeal et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	cachaca	110
2008	Cordero et al.	GC×GC-qMS	HS-SPME	Targeted profiling; Advanced fingerprinting	roasted coffee; roasted hazelnuts	8
2008	Klimankova et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	<i>Ocimum basilicum</i> L.	111
2009	Cajka et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	honey	92
2009	Cardeal & Marriott	GC×GC-TOFMS	HS-SPME	Targeted profiling	cachaca	112
2009	de Souza et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	cachaca	113
2009	Humston et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling; PARAFAC	cocoa beans	114
2009	Lojzova et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	potato chips	115
2009	Rochat et al.	GC-O; GC×GC-TOFMS	HS-SPME	Targeted profiling	shrimp aroma	78
2009	Vaz-Freire et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling; Advanced fingerprinting	EVO oil	65
2010	Breme et al.	GC-O; GC×GC-TOFMS	Direct injection	Targeted profiling	indian cress absolute	81
2010	Cordero et al.	GC×GC-qMS	HS-SPME	Targeted profiling; Advanced fingerprinting	juniper and roasted coffee	68
2010	Cordero et al.	GC×GC-qMS	HS-SPME	Targeted profiling; Advanced fingerprinting	<i>Corylus avellana</i>	70
2010	Humston et al.	GC×GC-TOFMS	HS-SPME	Untargeted profiling; PARAFAC	cocoa moisture damage	116
2010	Maikhunthod et al.	switcable GC×GC/MDGC-O	EO hydrodistillation	Targeted profiling	lavander eo	40
2010	Schmarr & Bernhardt	GC×GC-qMS	HS-SPME	Targeted profiling; Advanced fingerprinting	wine	66
2010	Schmarr et al.	GC×GC-qMS	HS-SPME	Targeted profiling	MOX red wine	67
2010	Schmarr et al.	GC×GC-qMS	SPE	Targeted profiling	wine	74
2010	Silva et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	marine salt	117
2010	Stanimirova et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	honey	118
2011	Chin et al.	GC-O/GC×GC-FID/FPD	SPE	Targeted profiling	wine, roasted coffee	119
2011	Gogus et al.	GC×GC-TOFMS	Direct TD	Targeted profiling	<i>Pistacia terebinthus</i>	60

2011	Pietra Torres et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	MLF wine	94
2011	Robinson et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	wine	84
2011	Robinson et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	wine	87
2011	Tranchida et al.	GC×GC-qMS; Capillary Flow Modulation	EO hydrodistillation	Targeted profiling	<i>Mentha spicata</i>	37
2011	Weldegergis et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	pinotage wines	120
2012	Chin et al.	GC-O/MDGC-FID/GC×GC-TOFMS	HS-SPME	Targeted profiling	shiraz wine	39
2012	Chin et al.	GC-O/GC×GC-FID; GC×GC-TOFMS	cumulative HS-SPME	Targeted profiling	shiraz wine	43
2012	Kiefl et al.	GC×GC-qMS	HS-SPME	Targeted profiling; Advanced fingerprinting	<i>Corylus avellana</i>	71
2012	Omar et al.	GC×GC-qMS/FID; Capillary Flow Modulation	FUSE	Targeted profiling	oregano; rosemary	121
2012	Villire et al.	GC-O/MS; GC×GC-TOFMS	Various HS approaches	Targeted profiling	cider	61
2012	Welke et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	merlot wine	77
2013	Bordiga et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling; Advanced fingerprinting	muscat wine	89
2013	Cordero et al.	GC×GC-qMS	Various HS and in-solution sampling	Targeted profiling; Advanced fingerprinting	dried milk	63
2013	Inui et al.	GC×GC-TOFMS	LLE	Targeted profiling	hop; beer	86
2013	Jelen et al.	GC×GC-TOFMS	SAFE	Targeted profiling	tempeh	80
2013	Kiefl et al.	GC×GC-TOFMS	SAFE	Targeted profiling	<i>Corylus avellana</i>	64
2013	Kiefl & Schieberle	GC×GC-TOFMS	SAFE	Targeted profiling	<i>Corylus avellana</i>	85
2013	Langos et al.	GC×GC-TOFMS	SAFE	Targeted profiling	beer	122
2013	Maikhunthod & Marriott	GC-O/GC×GC-FID; GC×GC-TOFMS	HS-SPME	Targeted profiling	dried spice	44
2013	Majcher et al.	GC×GC-TOFMS	SAFE	Targeted profiling	cereal coffee	79
2013	Mommers et al.	tunable-GC×GC-TOFMS	HS-SPME	Targeted profiling	roasted coffee	45
2013	Nicolotti et al.	GC×GC-qMS	MHE-SPME	Targeted profiling; Quantitative fingerprinting	<i>Corylus avellana</i>	59
2013	Rivellino et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	honey	123
2013	Samykanno et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	strawberry	88
2013	Tranchida et al.	LC-GC×GC-qMS	EO cold pressing	Targeted profiling	citrus essential oil	38
2013	Van Der Wat et al.	GC×GC-TOFMS; GC-O	Multi-channel PDMS traps	Targeted profiling	rosemary	62
2013	Willner et al.	GC×GC-TOFMS	SAFE	Targeted profiling	brandy	124
2013	Zhang et al.	GC×GC-TOFMS	SDE	Targeted profiling	Tea	96
2014	Bernal et al.	GC×GC-FID Capillary Flow Modulation	S-HS	Advanced fingerprinting	roasted almonds	36
2014	Bordiga et al.	GC×GC-TOFMS	HS-SPME	Targeted profiling	wine	125
2014	Dugo et al.	GC×GC-qMS/FID	HS-SPME	Targeted profiling	wine	126
2014	Purcaro et al.	GC×GC-qMS	HS-SPME	Targeted profiling; Advanced fingerprinting	EVO oil	72

Figure 1

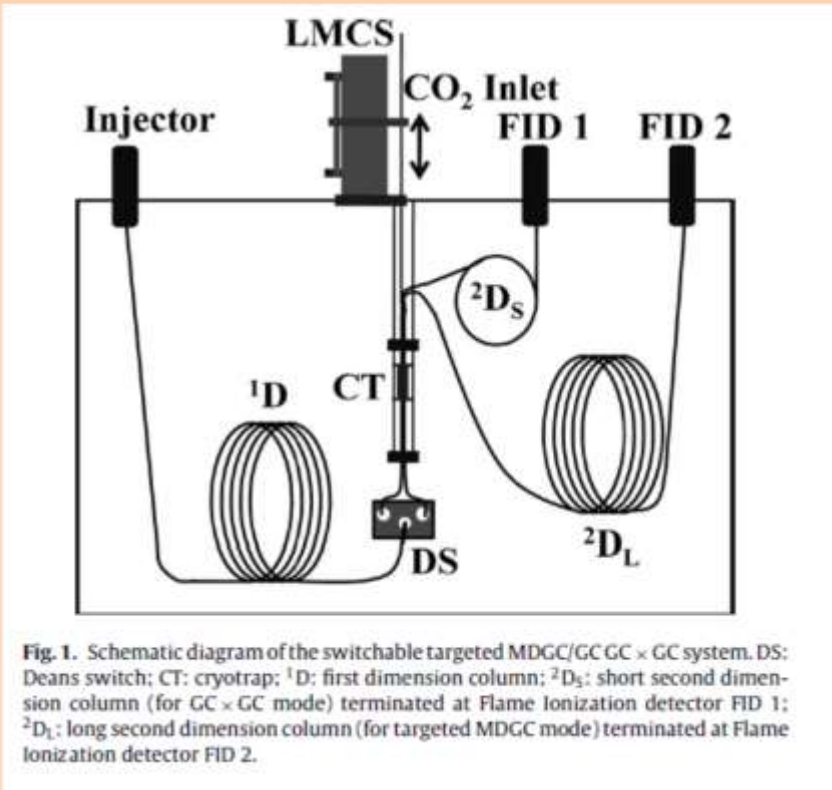

1205

1206

Figure 2


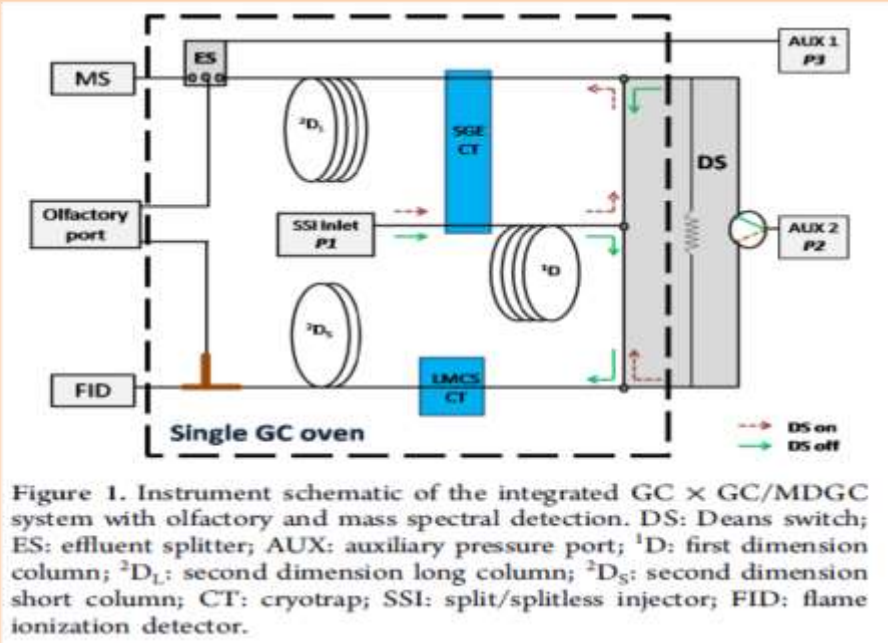



Figure 3

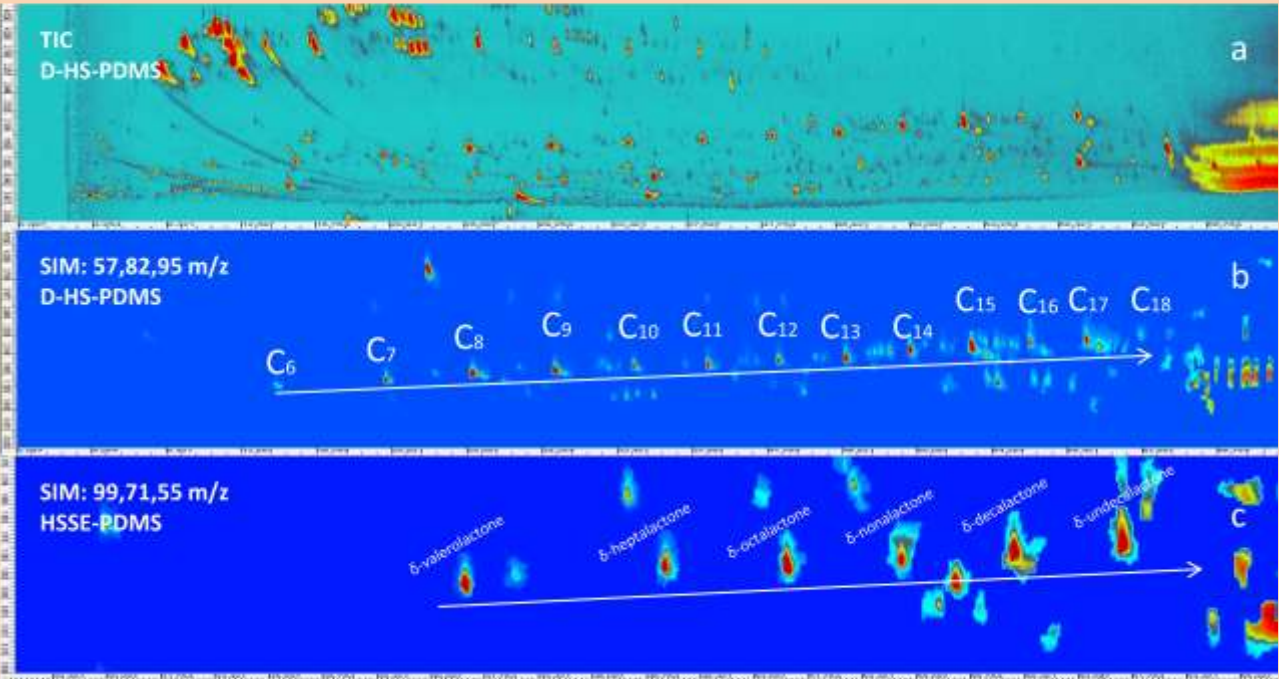

1209

1210


Figure 4

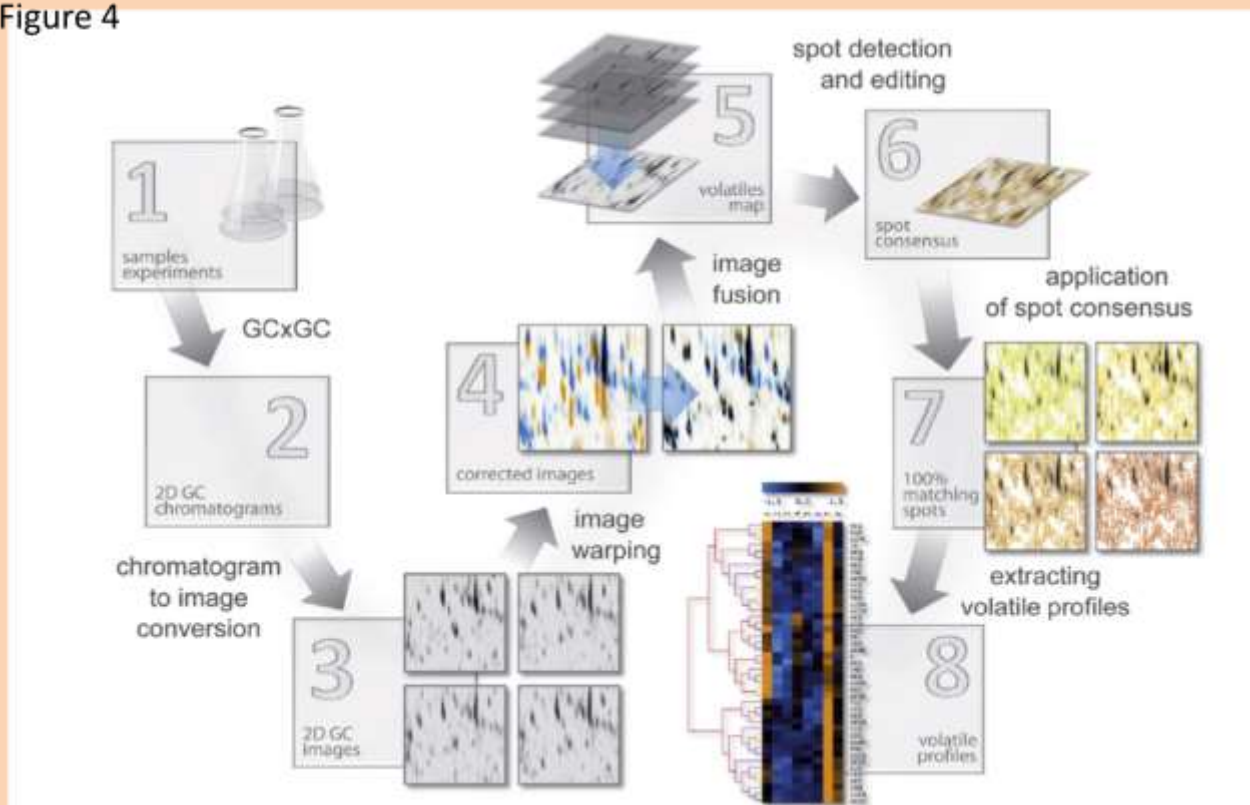

1211

1212

Figure 5

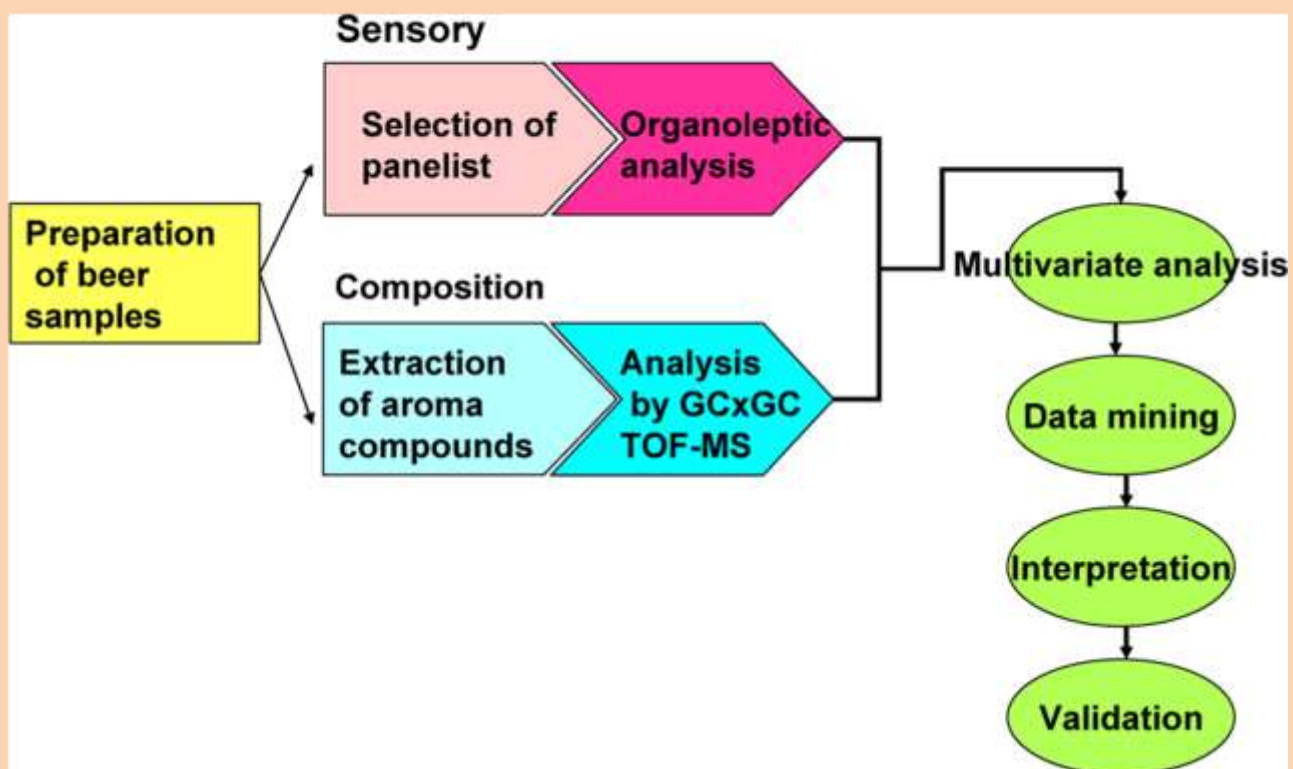

1213

1214

Figure 6

